

Rec'd PCT/PTO 18 JAN 2003

18/500988 IL 03/00026

REC'D 22 MAY 2003

WIPO

PCT

PA 996918

# THE UNITED STATES OF AMERICA

**TO ALL TO WHOM THESE PRESENTS SHALL COME:**

**UNITED STATES DEPARTMENT OF COMMERCE**

**United States Patent and Trademark Office**

**April 18, 2003**

**THIS IS TO CERTIFY THAT ANNEXED HERETO IS A TRUE COPY FROM THE RECORDS OF THE UNITED STATES PATENT AND TRADEMARK OFFICE OF THOSE PAPERS OF THE BELOW IDENTIFIED PATENT APPLICATION THAT MET THE REQUIREMENTS TO BE GRANTED A FILING DATE UNDER 35 USC 111.**

**APPLICATION NUMBER: 60/391,575**

**FILING DATE: June 27, 2002**

**BEST AVAILABLE COPY**

**PRIORITY DOCUMENT  
SUBMITTED OR TRANSMITTED IN  
COMPLIANCE WITH  
RULE 17.1(a) OR (b)**

**By Authority of the  
COMMISSIONER OF PATENTS AND TRADEMARKS**



*Gloria A. Murray*  
**GLORIA A MURRAY**  
**Certifying Officer**

27/06 '02 14:12

0872 3 5825554

DR. M. FRIEDMAN \*\*\* BILL

@008/011

FROM : IMT LTD

FAX NO. : 4972 8 9409988

JUN 27 2002 02:19 PM

*WJD*

### ASSIGNMENT

For good and valuable consideration, the receipt and sufficiency of which is hereby acknowledged, the undersigned:

Amir Arav

Meir Uri

(hereinafter called the "assignor(s)"), hereby sell(s), assign(s) and transfer(s) to:

Interface Multigrad Technology Ltd.  
3 Hamazmera St.  
Ness Ziona 70400  
Israel

(hereinafter called the "assignee(s)"), its/his successors, assignees, nominees or other legal representatives, the Assignor's entire right, title and interest in and to the invention entitled:

### **EMBRYO GUARD**

described and claimed in the following patent applications:

U.S. Provisional Application identified as Attorney docket No. 791/14 and executed the same date as this assignment;

and in and to said Patent Applications, and all original and reissued Patents granted therefor, and all divisions and continuations thereof, including the right to apply and obtain Patents in all other countries, the priority rights under International Conventions, and the Letters Patent which may be granted thereon;

Signed and sealed this 27 day of June 2002

*[Signature]*  
Amir Arav

*[Signature]*  
Meir Uri

UNITED STATES PATENT AND TRADEMARK OFFICE  
DOCUMENT CLASSIFICATION BARCODE SHEET




# As-Filed New Application

Level -1  
Version 1.1  
Updated - 8/01/01

UNITED STATES PATENT AND TRADEMARK OFFICE  
DOCUMENT CLASSIFICATION BARCODE SHEET



Transmittal

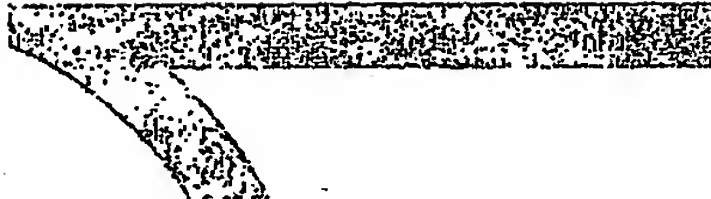

Level - 2  
Version 1.1  
Updated - 8/01/01

60391575.062702

A/PROV

Please type a plus sign (+) inside this box

# PROVISIONAL APPLICATION FOR PATENT COVER SHEET

This is a request for filing a PROVISIONAL APPLICATION FOR PATENT under 37 CFR 1.53 (c).

 11129 U.S. PTO  
 06/27/02

 60/391575  
 U.S. PTO  
 06/27/02

DocId Number: 741/14

INVENTOR(S)/APPLICANT(S)		Residence (City and State or Foreign Country)	
Class Name (last and middle if any)	Family Name or Surname		
AMIE MEIR	ARAY URI	BET DAGAN, ISRAEL	
Additional Inventors are being named on page 2 attached hereto		BET HASHITA, ISRAEL	
TITLE OF THE INVENTION (250 characters max)			
EMBRYO GUARD			
Direct all correspondence to:		CORRESPONDENCE ADDRESS	
<input type="checkbox"/> Customer Number			
<input type="checkbox"/> OR		Place Customer Number Bar Code Label Here	
Firm or Individual Name	Mark M. Friedman		
Address	9003 FLORIN WAY		
City	UPPER MARIBORO	State	MD
Country	US	Telephone	301-9521011
		ZIP	20732
		Fax	301-9529003
ENCLOSED APPLICATION PARTS (check all that apply)			
<input checked="" type="checkbox"/> Specification	Number of Pages	Applicant is Small Entity	
<input type="checkbox"/> Drawing(s)	Number of Sheets	<input checked="" type="checkbox"/> Other (specify) Assignment	
METHOD OF PAYMENT OF FILING FEES FOR THIS PROVISIONAL APPLICATION FOR PATENT (check one)			
<input type="checkbox"/> A check or money order is enclosed to cover the filing fees			
<input checked="" type="checkbox"/> The Commissioner is hereby authorized to charge filing fees or credit any overpayment to Deposit Account Number:			
The invention was made by an agency of the United States Government or under a contract with an agency of the United States Government.			
<input checked="" type="checkbox"/> No.			
<input type="checkbox"/> Yes, the name of the U.S. Government agency and the Government contract number are:			
Respectfully submitted,			

 SIGNATURE Mark M. Friedman DATE 27 Jun 02  
 TYPED or PRINTED NAME Mark M. Friedman REGISTRATION NO. 33,883  
 (If appropriate)  
 TELEPHONE (703) 415-1581


 USE ONLY FOR FILING A PROVISIONAL APPLICATION FOR PATENT  
 SEND TO: Box Provisional Application, Assistant Commissioner for Patents, Washington, DC 20231

(Page 1 of 1)

UNITED STATES PATENT AND TRADEMARK OFFICE  
DOCUMENT CLASSIFICATION BARCODE SHEET



# Specification

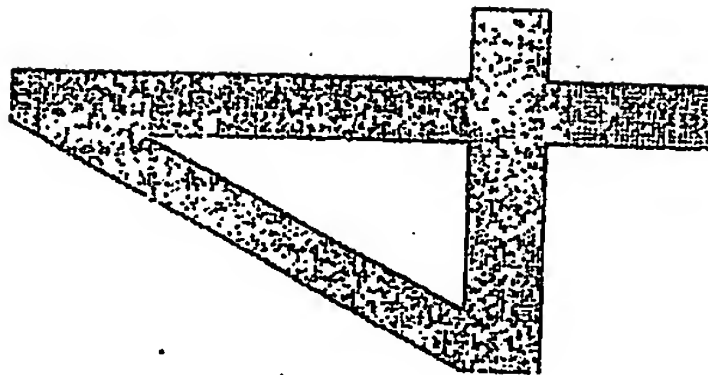

Level 2  
Version 1.1  
Updated - 8/01/01



June 2002

IMT Ltd.

## Microscopic Monitoring

**Confidentiality Note**  
This document contains proprietary  
information belonging to IMT Ltd.  
Reproduction of this document, disclosure of  
its content or any other use of information  
herein is forbidden unless specifically  
authorized by IMT Ltd. in writing.

## Technology and Products



## Microscopic Monitoring

During its many years of research in the fields of cryopreservation and reproduction, IMT has developed several techniques and supporting technologies to improve its research process. The MMS technology is a result of one of these developments.

Initially, the company's researchers needed a tool that allows microscopic monitoring of cells during the freezing process. The solution was a miniature device, consisting of a CCD camera, a microscopic long-distance objective lens and a special adapter. That was the beginning of a new field for IMT and based on that simple device, the company develops several important products, including the EmbryoGuard.

### 1. The Technology

Microscopic monitoring is very common in biology, as it is a basic tool for most procedures. In every lab there are usually several types of microscopes and most of the devices have a special adapter for video or CCD camera.

However, this simple approach has several limitations. First, the image must go through the optical system of the microscope. This means that the microscope should be a high quality device, so that the quality of imaging remains untouched. Such high quality microscopes are very costly. In addition, the size of the microscope limits its applications. If one wants to monitor cell development inside an incubator, he can "build" an incubator that surrounds an existing microscope, or he can take the cells out of the incubator and place them under a microscope in a warm environment.

IMT has succeeded in minimizing the whole system, while avoiding the need for a microscope. A very small CCD or video camera, with a special adapter and a microscopic lens, can provide the same results as an expensive and high quality system.

An improvement to this basic technology was introduced when the Company developed a robotic system that can control multi-sample monitoring, with X-Y-Z micro movement. This system is based on one or two microscopic CCD cameras, together with robotic features that can move the cameras and the samples. This system is then placed inside a standard incubator. A control unit placed outside the incubator helps to control focusing and illumination, changing between samples, and more. The



control unit can be handled manually or by computer software and the images are screened on a standard monitor.

## 2. Products and Applications

The Company is developing several products based on that technology. The main one is the Embryo Guard.

The Embryo Guard (EG) is a robotic system for microscopic monitoring and control over embryo development during IVF procedure. It has computerized control and software that assists with embryo evaluation, supports the selection process, and controls the matching process (patient/oocyte/sperm/embryo).

During IVF procedure, it is extremely important to monitor the development of fertilized eggs from the moment of fertilization up to the stage when 2-3 embryos are selected for transplantation. The importance of this monitoring is derived from the fact that eventually, the best embryos should be chosen for transplantation. Therefore, the clinician must watch very carefully every development stage of the embryos.

### Problems being addressed

Change in conditions - Currently, the monitoring procedure is done manually, by taking the embryos out of the incubator, placing them under a microscope and investigating their development.

This approach has several disadvantages that usually damage the embryos. To assure the best conditions for embryo development, it is essential that the embryo remains in a stable controlled environment, as provided by the incubator. Any change in these conditions can easily harm the embryo. Therefore, the procedure of taking the embryo out of the incubator, although it is necessary, has a bad effect on the embryo development. In addition, the optimal way to evaluate embryo development is to monitor it every 3 hours, but again, since this monitoring might be too risky, most IVF labs prefer to perform this evaluation much less frequently. Another problem raised from the need to monitor an embryo under a microscope is that the embryo must be under a special solution (oil) that can damage it.

The EG solves these problems by providing continuous monitoring of the embryo, without taking it out of the incubator. The EG, automatically, monitors each embryo

every 3 hours, or continuously (time-laps recording) and stores this data (as image files) on the embryo records. This process is done inside the incubator which means that embryos do not experience any change in condition and also, there is no need to use oil or any other solution that may cause damage.

Another important advantage of the EG is the ability to control it from a distance, using the Internet. The embryo specialist does not have to be present in the IVF lab, and he can control the whole process using a standard computer connected to the Internet.

Standards and data recording - Another problem addressed is the lack of standards and data recording. Currently, there is no software application that supports the IVF procedure in terms of embryo development and selection of eggs and embryo. Without such supporting software, IVF labs collect data on the embryos in a variety of ways, with no specific standard or quality control, and in addition, most of this data collection is done by paper work.

The EG includes a software that satisfies these needs. Each embryo in the incubator has its own record, containing all the information from the initial stage. The software also automatically collects and stores pictures of each embryo in each stage. In addition to data collecting, the software also helps to evaluate the embryo by indicating in which stage it should be, how many cells it should have, what should be the next stage and the timing for this stage, etc. It can also provide a multi-embryo screen that helps to compare their visual shape.

Matching - IVF labs pay a lot of attention to the issue of matching between oocytes and sperm, or embryos and patient. Even a minor mistake could be a personal disaster for the future parents and a major legal problem to the lab. One component of the EG aims to solve this problem. The EG has a unique matching system that makes sure that no mistake can happen. Every sperm sample and every oocyte are labeled (on the dish) with a barcode labeling system. This procedure is then stored in the computer record as the first step of the IVF procedure. From the moment of labeling, every procedure must first go through a barcode reader that stores the information under the patient record. Before fertilization, the EG software identifies both patients and indicates if there is matching or not. The fertilized eggs are then placed inside the incubator (again - after going through the barcode reader). If during the incubation

period, or even before transplanting the embryos, the clinician needs to take the embryos outside the incubator, the software identifies the specific dish and let the clinician take only this one.

#### Schedule

IMT plan to introduce the first commercial version of the EG on July 2002.

#### Consumable products

The Embryo Guard can handle up to 12(?) dishes simultaneously. Each dish is for one specific patient and it can contain up to 10(?) embryos. The dish is a regular and standard dish, sterilized, and is currently available in the market.

For locating and identifying the embryos in specific locations around the dish, the company offers a special sticker for each dish, which also contains barcode ID. The EG cannot operate without IMT's stickers.

embryo guard notes

## Embryo Guard

Notes from interview with Amir Arav, 26 June 27-06-2002

Based on US 6,166,761. Now full robotic system.

Point 1: Follow development of embryos because in vitro fertilization depends on first cleavage. Need to know timing of cleavage. Implant three-day-old embryos. Hard to tell which one to pick.

### Point 2:

Prior art: embryos removed from incubator to exchange medium.  
New: automatic medium exchange. Add medium, remove medium, or do both. Gas

the medium before or after warming up the medium.

### Point 3:

Zona pellucida stays thick in an incubator. As a result, the embryo may die. Focus a laser beam through a microscope to heat the zona pellucida to cut it open. Do this inside the incubator. Cutting may be done manually by technician, or automatically. Prior art is to remove embryo from incubator and cut zona pellucida outside incubator.

Point 4: Fluorescent markers. Also for preimplantation diagnostics.

Point 5: Insemination inside incubator.

Point 6: Preparation for cryopreservation: is a special case of Point 1.

In general: micromanipulation of oocytes and embryos is done inside the incubator.

Matching using bar code (or equivalent: remotely readable chip, imaging, etc.)  
Identify gametes, oocytes, sperm upon collection.

Stickers on test tubes, vials, petri dishes (containers generally). ID text (parents names etc.) matched automatically to ID code.

Box outside incubator has place for one test tube and one petri dish.  
Match test tube to petri dish based on bar code. Match embryos and cryogenic vials

(need liquid-nitrogen-resistant bar code) for cryopreservation. Match again when transferring embryos to womb.

Gamete Intra Fallopian Transfer. Zygote Intra Fallopian Transfer. Match at Pre-Implantation Diagnosis.

Management software: tells you what to do when (timing is critical). Warns if embryos are outside incubator too long. Collect history automatically.

embryo guard notes

To be able to handle 12 petri dishes in the same incubator:

Automatic orientation

CCD camera goes from drop to drop automatically

Search for embryo automatically or go to center of drop

Embryos can move, so image at least 10 times the area of an embryo

Digital magnification. Use high resolution CCD.

The Embryo Guard also can be used for other applications which require culturing cells or tissues in an incubator for a long period of time and to monitor the cells or tissue without removing them from the incubator.

Patent on EmbryoGuard

1. On line monitoring, time lags reordering, medium change over and assisted hatching of embryos inside the incubator.

It is been well recognized that the timing of the first cleavage and the morphology of the embryos determine the successful of the IVF procedure in addition, change over of medium and assisted hatching are other reasons for removal the embryos from the incubator which could be optimized if they could be performed inside the incubator. Many opening of the incubator i.e. for microscopic evaluation, medium change-over and assisted hatching, will affect the embryos culture condition (temperature, gas concentration and humidity).

We describe here a robotic system which will operate in the incubator with the following feature:

1. A microscopic follow up of the development of embryos inside the incubator with the following possibility: A. A real-time evaluation of the embryos using up to 4 different microscope CCDs which could be operate on 3D movements also by using internet compatibility. B. automatic photographing system for image of up to 12 different dishes in which 12 drops are placed in each of the Petri dish.
2. Medium change over is done by computerized injection to each of the drops with small volume of 1-10 microliter of fresh medium which is maintained cooled before injection and then warmed up, gassed and added to the drop.
3. Assisted hatching is done by laser beam performing on the zona pelucida in order to assist blastocyst hatching.
4. Evaluation using fluorescent markers which are loaded with the injectors and are detected by the embryoGuard.
5. Insemination could be done stepwise by the injectors insert sperms directly to the drops.
6. preparation oocytes or embryos for cryopreservation inside the drops using the computerized injector in a stepwise manner and according to the osmotic behavior of the cells.

## 2. Control matching using barcode system.

It is estimate that there are hundreds of mistakes in IVF matching worldwide. Identification of embryos is done by the technicians created humanity mistakes. We propose of using a adhesive sticker with a barcode for test tube and Petri dish. The EmbryoGuard read the barcode and identify the oocyte for matching with sperm. In a separate apparatus which is placed outside of the incubator (an EmbryoSeater). The matching could be performed in several levels:


1. inside the incubator before the oocytes are fertilized.
2. when sperm arrive to the lab
3. between sperm and oocytes
4. for PGD
5. for cryopreservation

In a case of no matching the EmbryoGuard will not aloud to be open and remove out the oocytes or embryos or done any other function.









2/08 '02 14:14 6872 3 682554 DR. B. FRIEDMAN +-+ BILL 2008/011

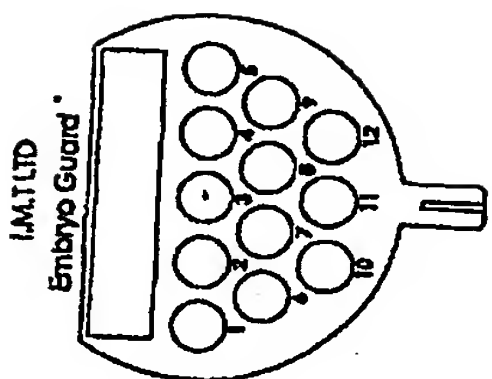


# The EMBRYOGUARD

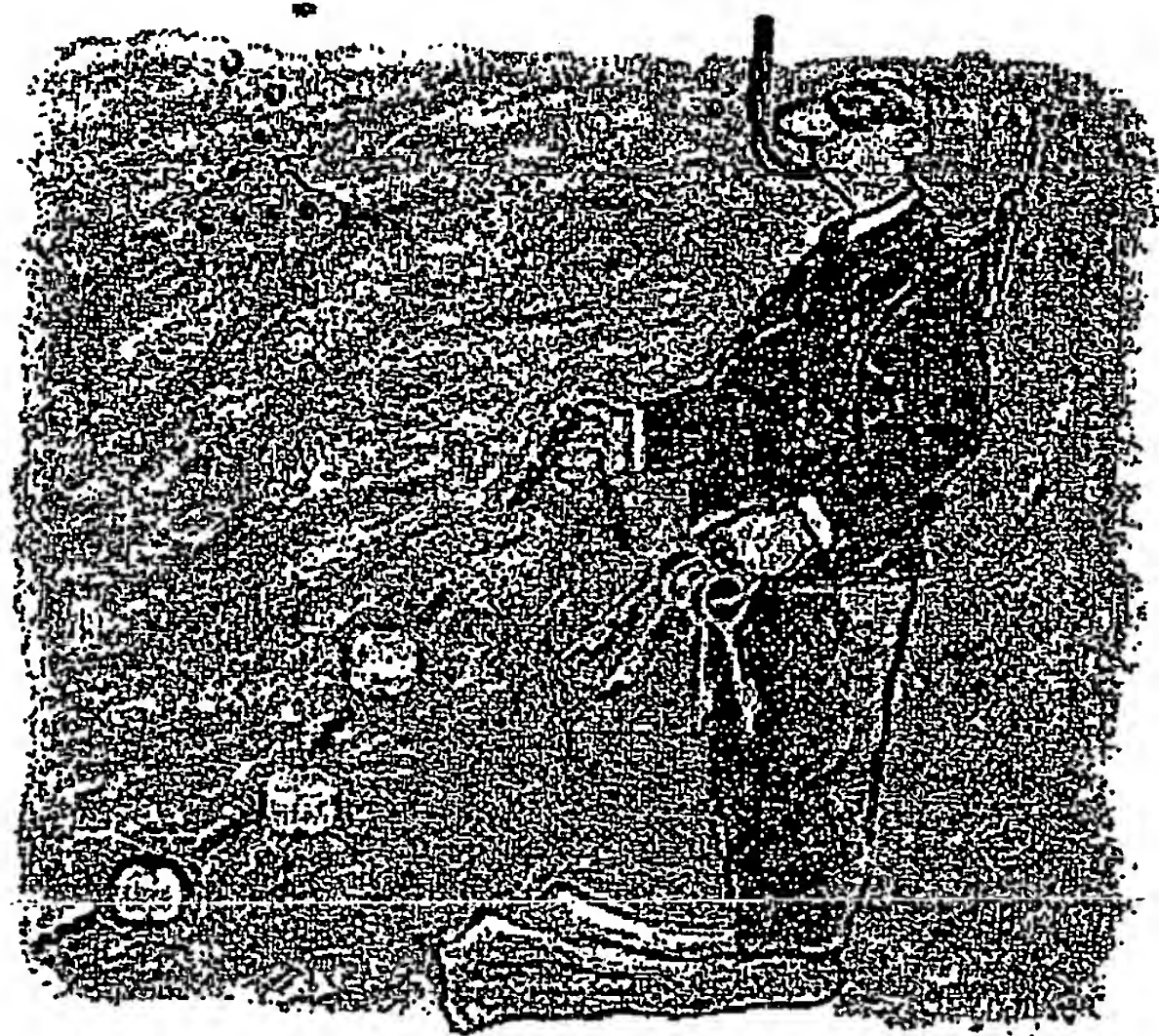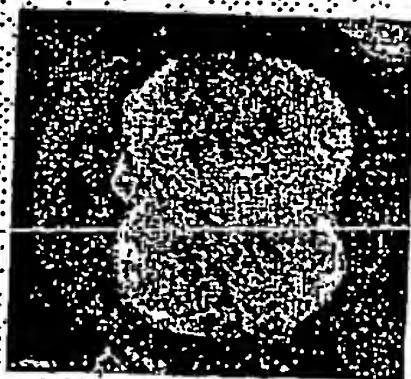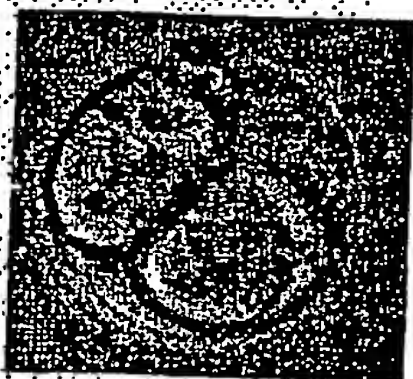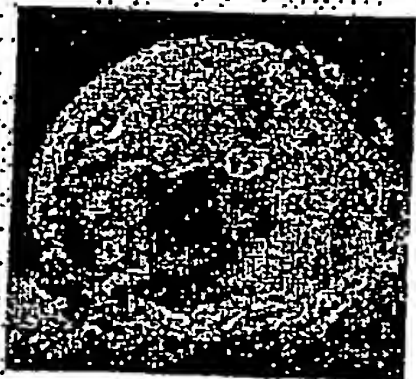


## The Embryoguard

- A new system which contains three microscopic CCD cameras located inside the IVF incubator, including active matching management tool

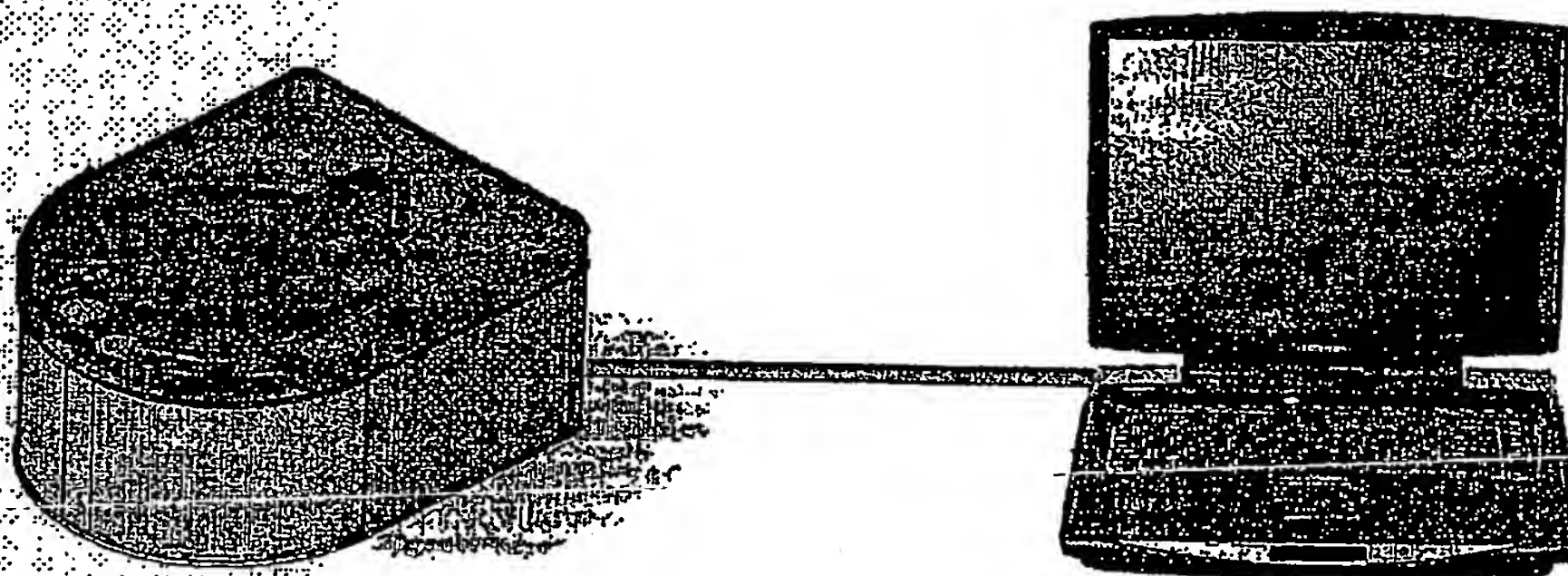





# **On-line monitoring & time-laps evaluation of embryos inside incubator**

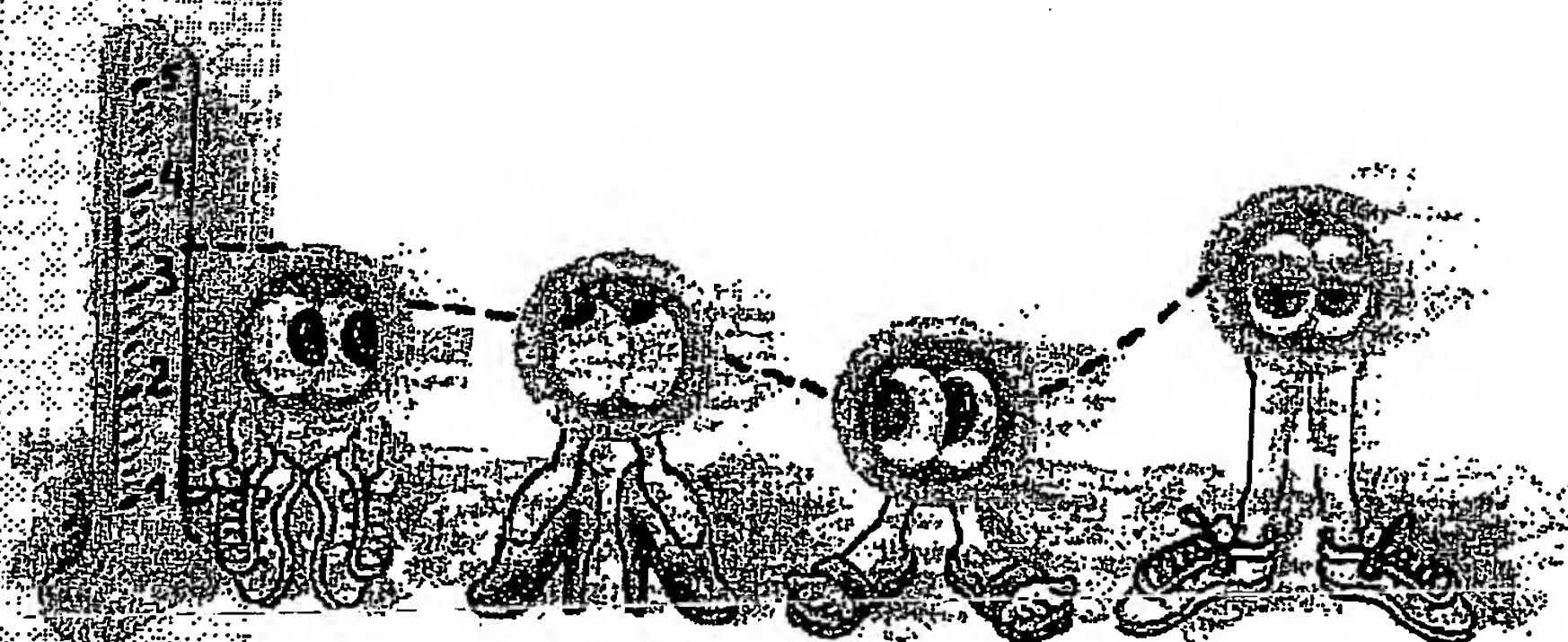

# **Improved success rates by selecting of embryos based on cleavage timing**

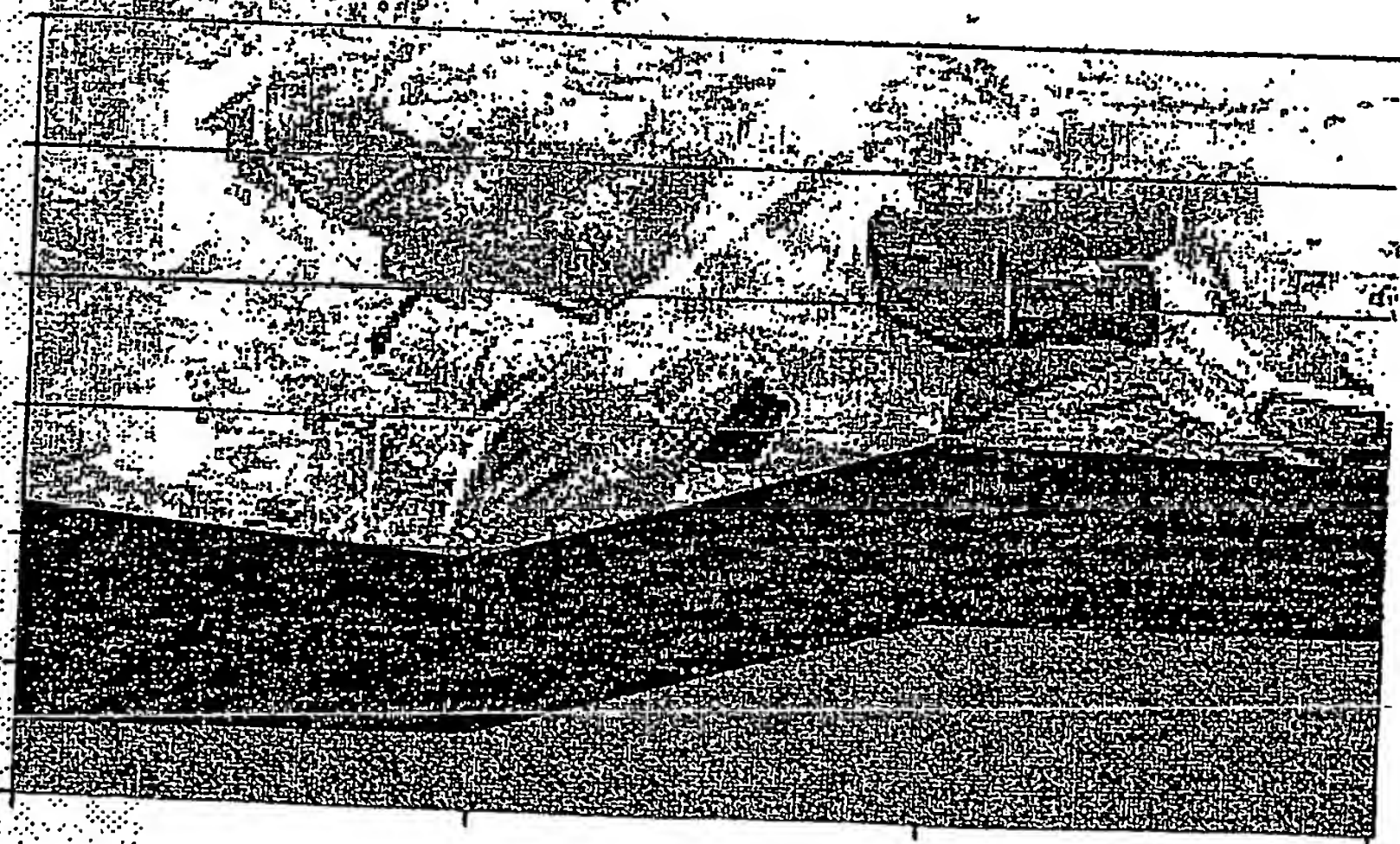


# **Optimal management of IVF lab procedures**

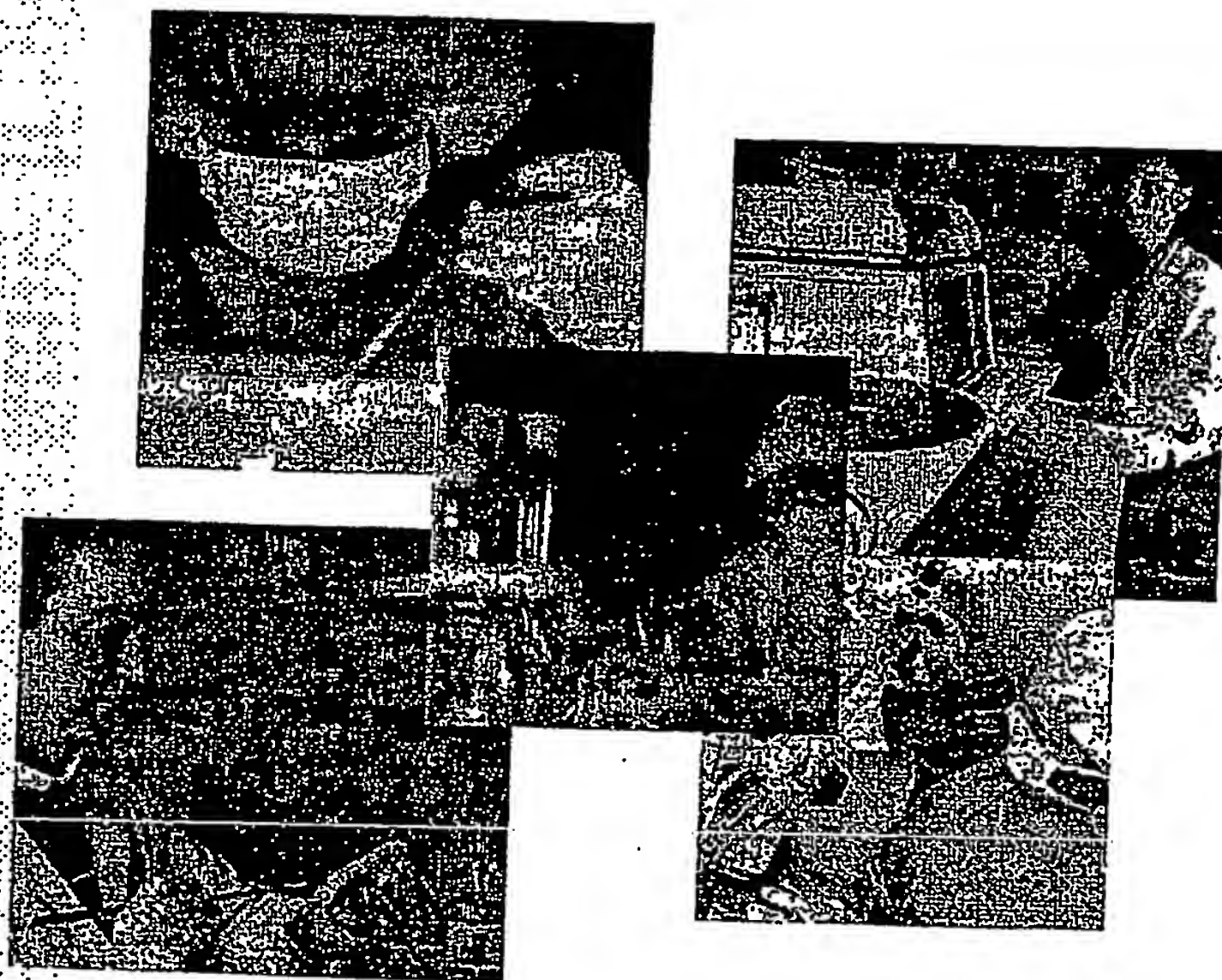





## Controlled matching utilizing the Barcode identification system

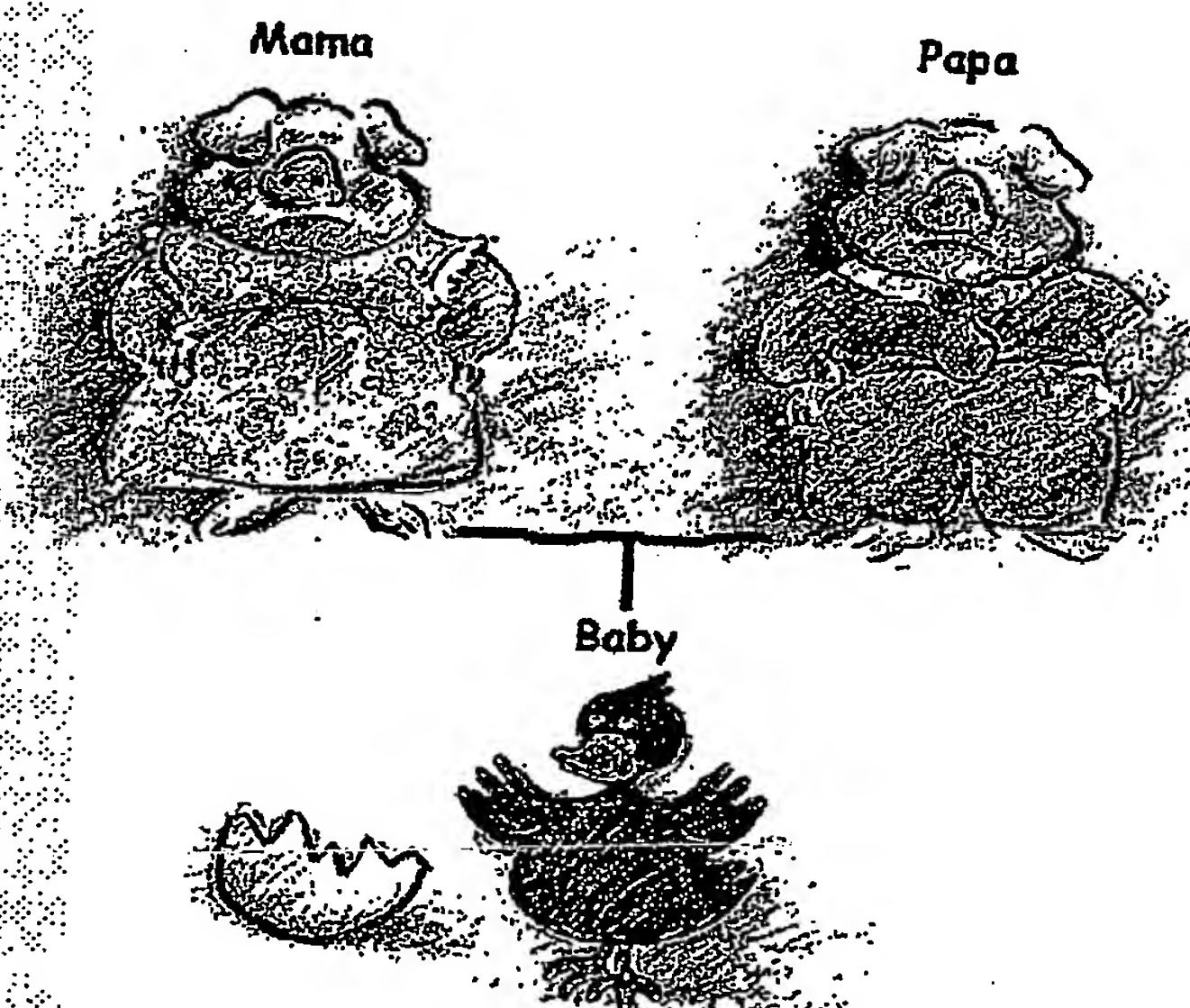

# Complete documentation

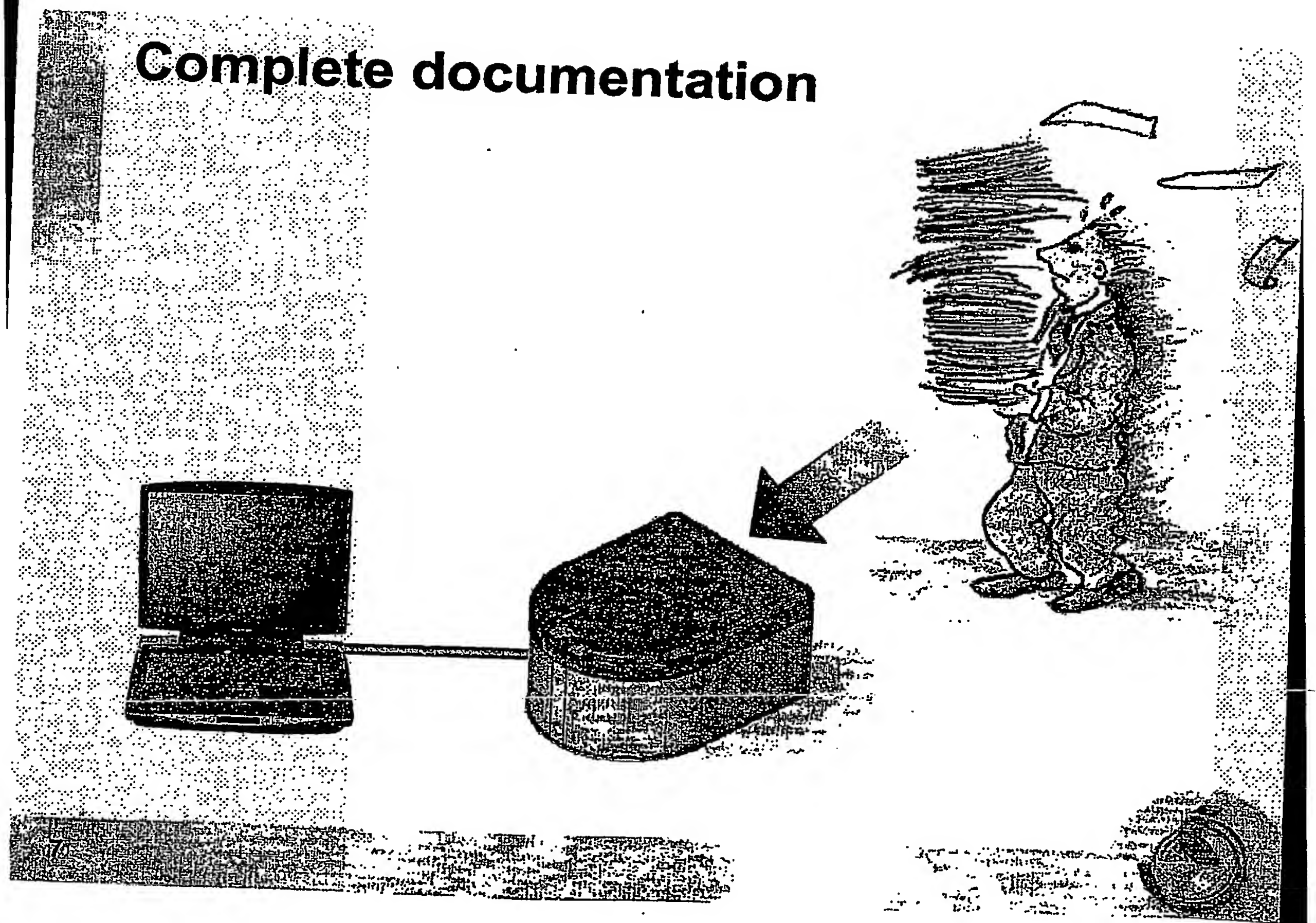


## On line control per procedure



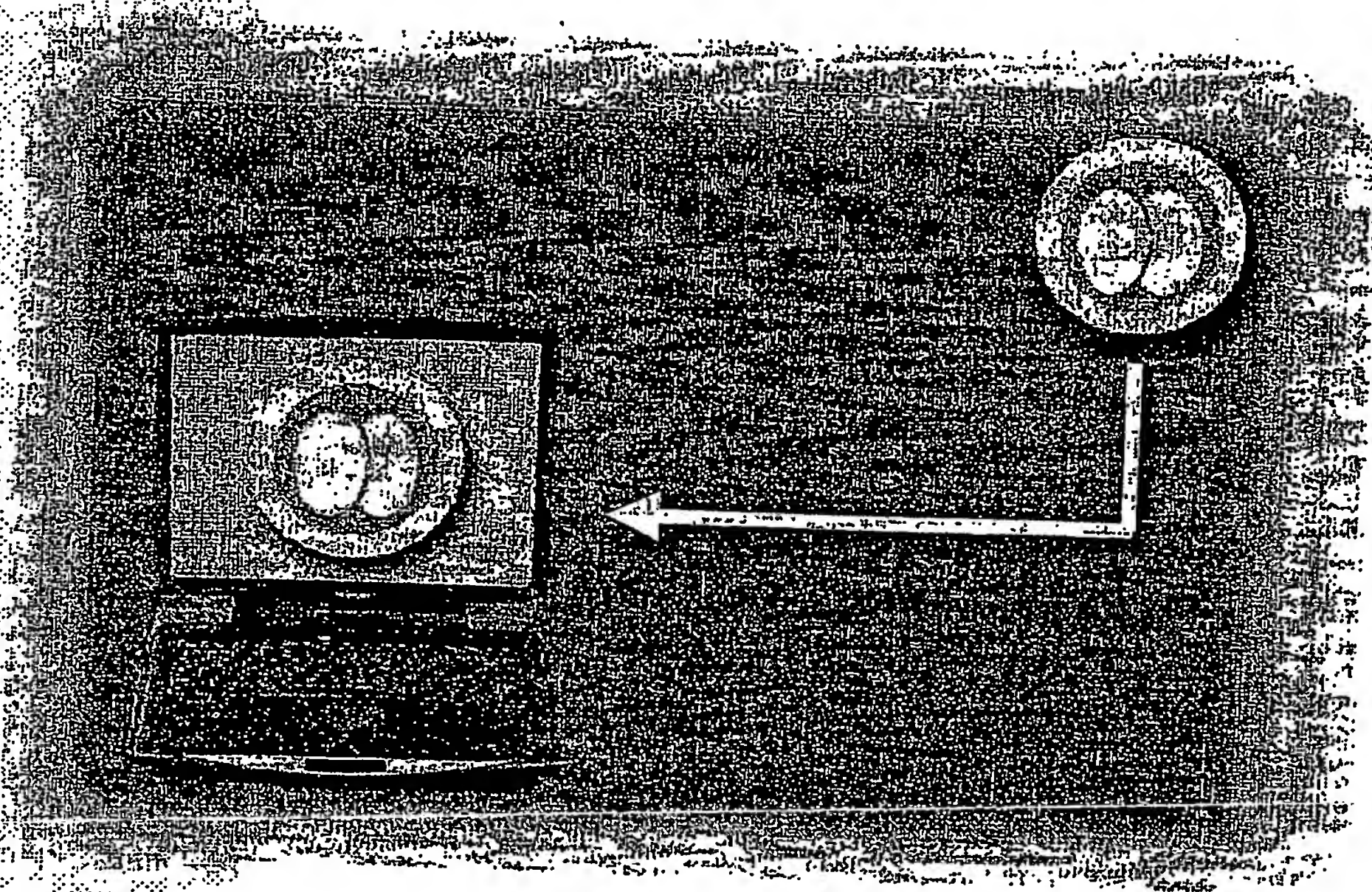




## Internet compatibility

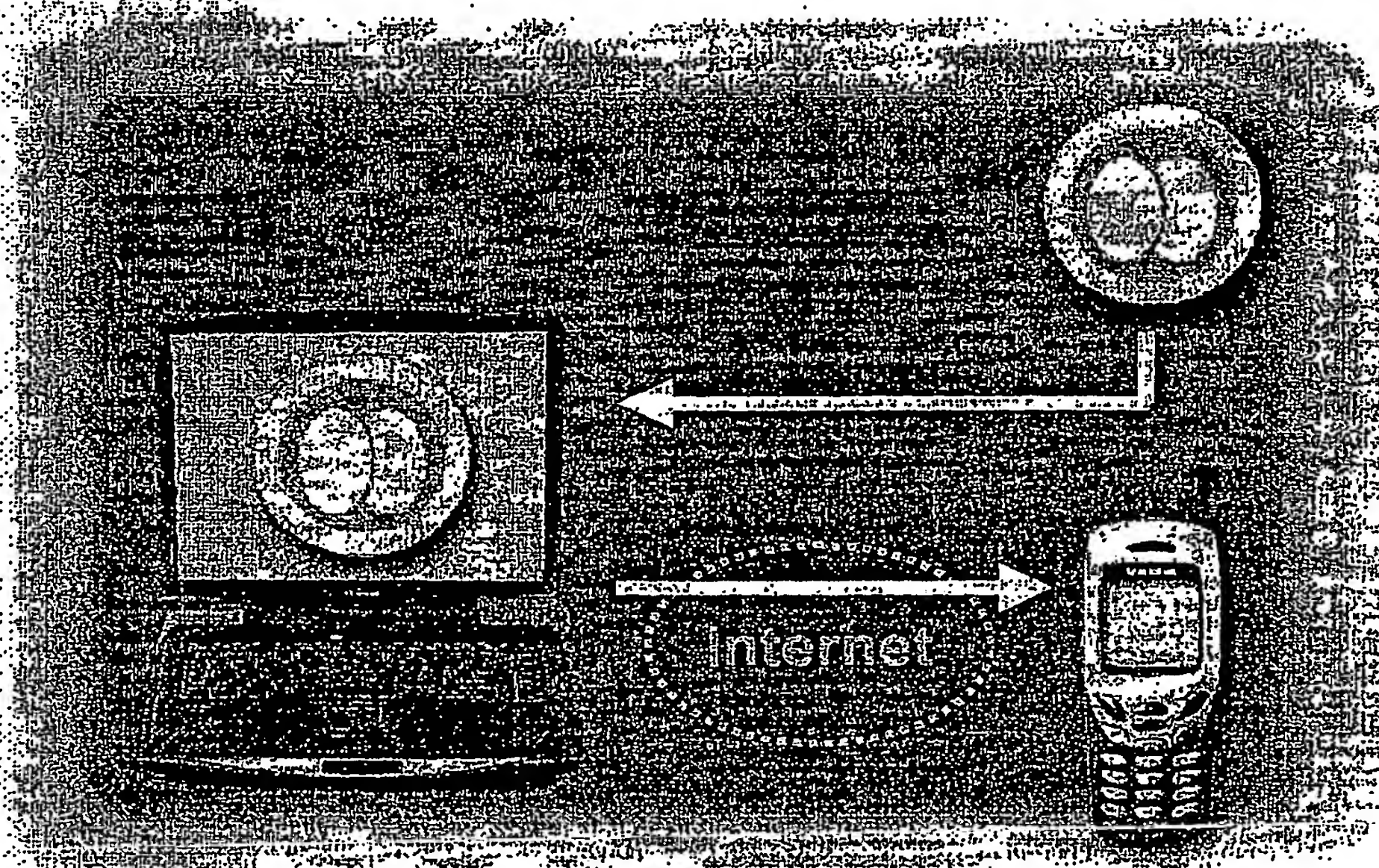


# Real time evaluation of up to 12 dishes simultaneously

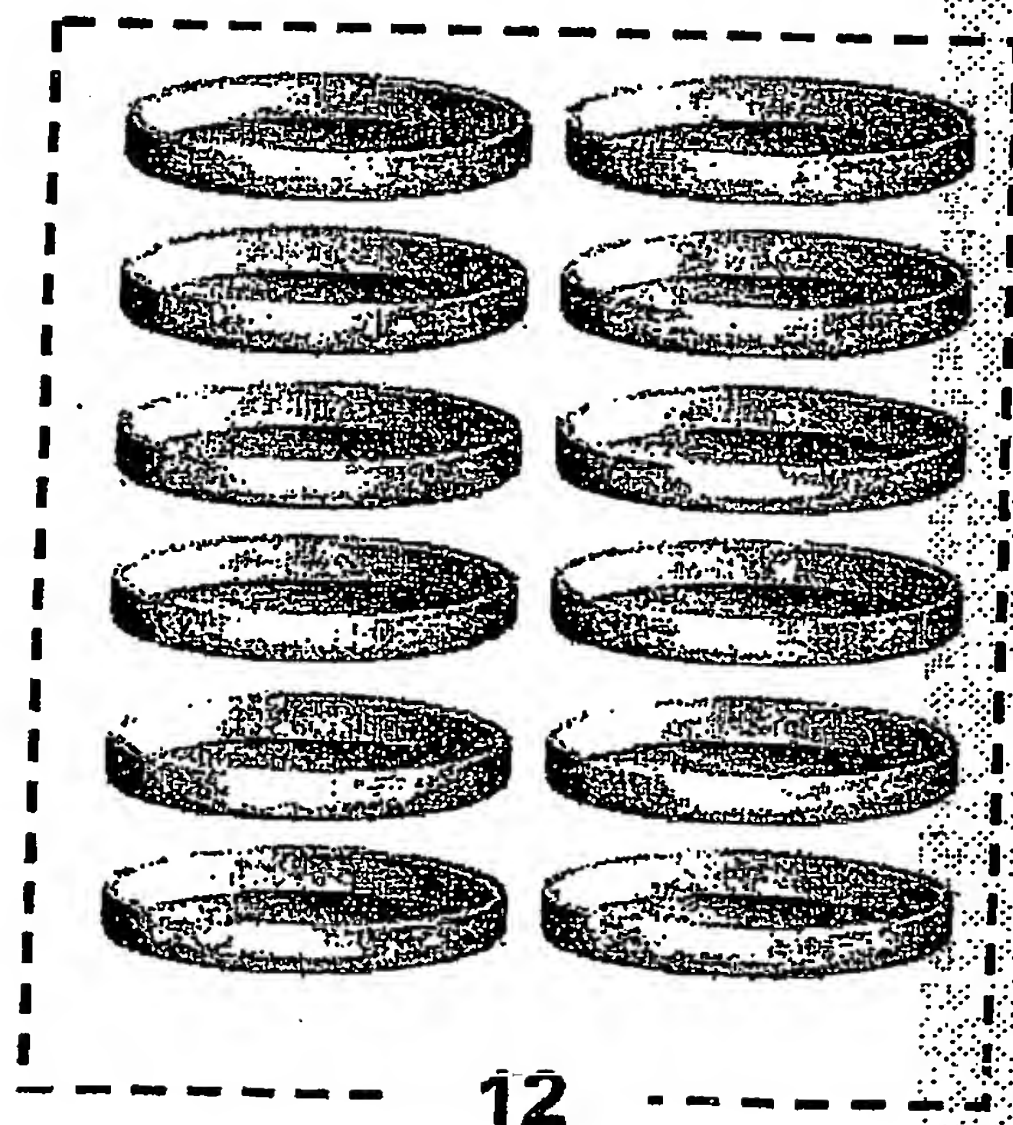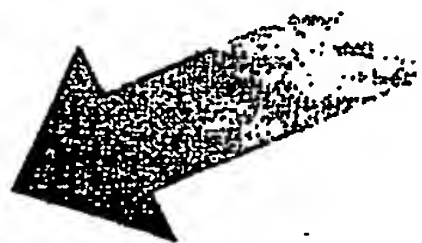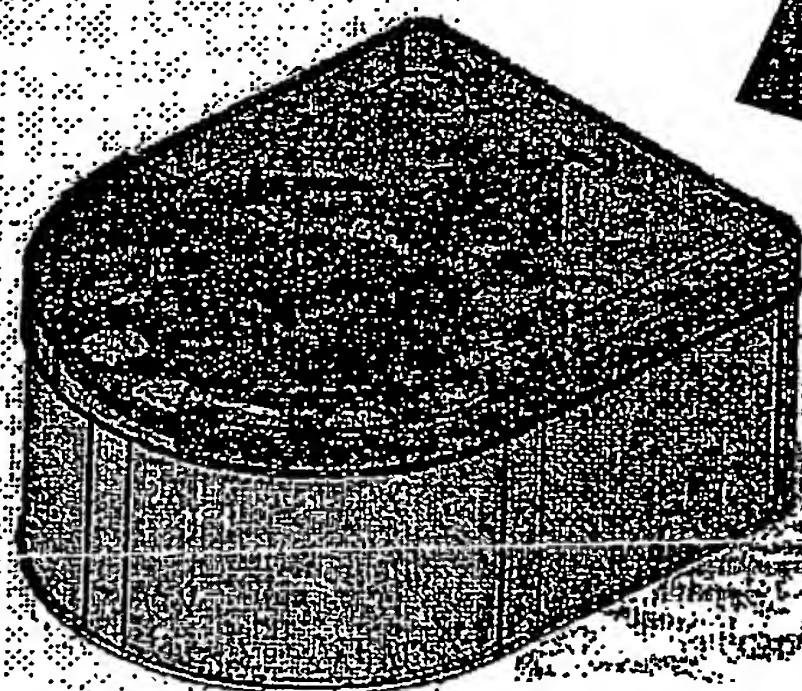

12




# **Full robotic 3D movement of CCD microscopic cameras inside the incubator**

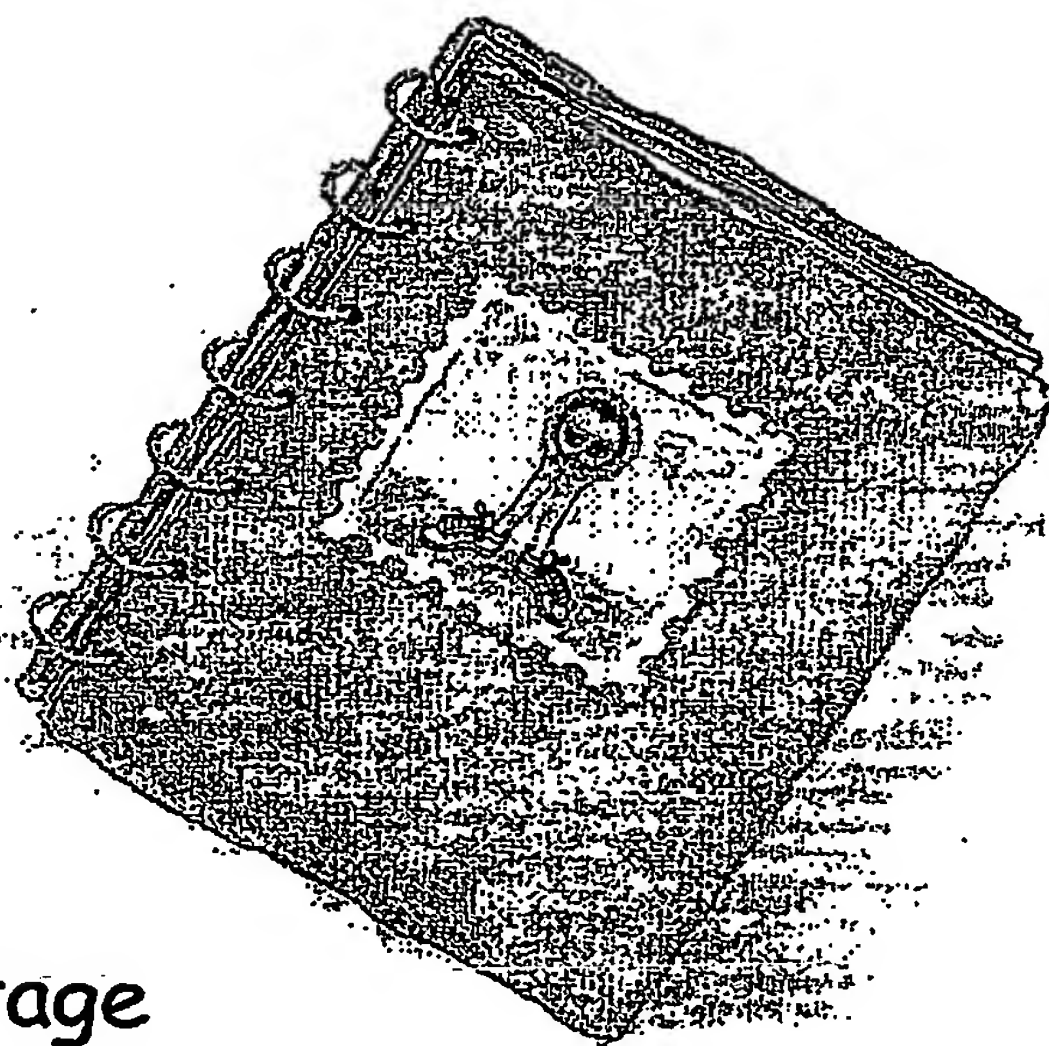

**Start your  
photo album  
from your 2PN stage**


**Reduces the need opening the incubator**

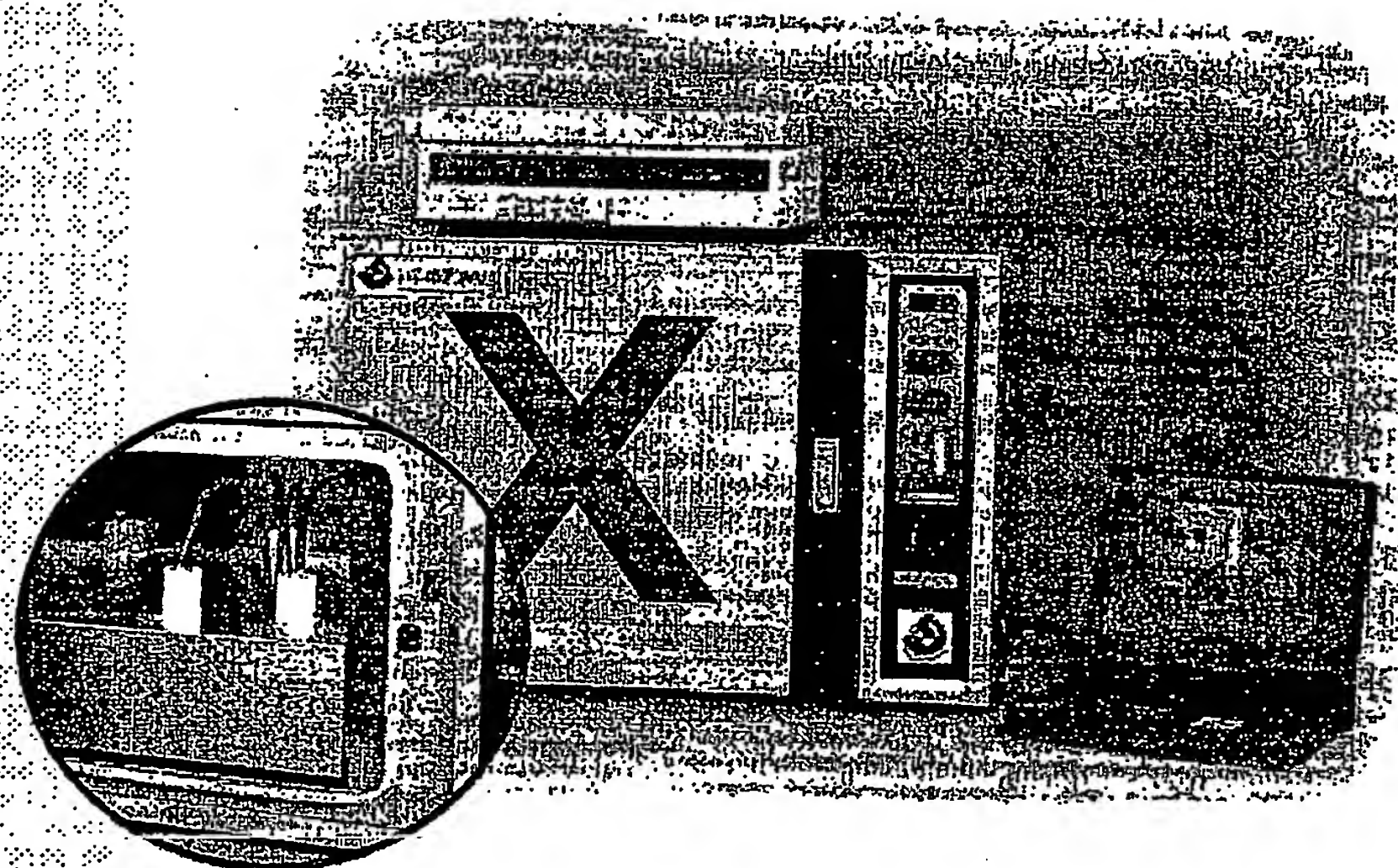





## For More Information


- E-mail: [embryoguard@cryo-imt.com](mailto:embryoguard@cryo-imt.com)
- Site: [www.cryo-imt.com](http://www.cryo-imt.com)

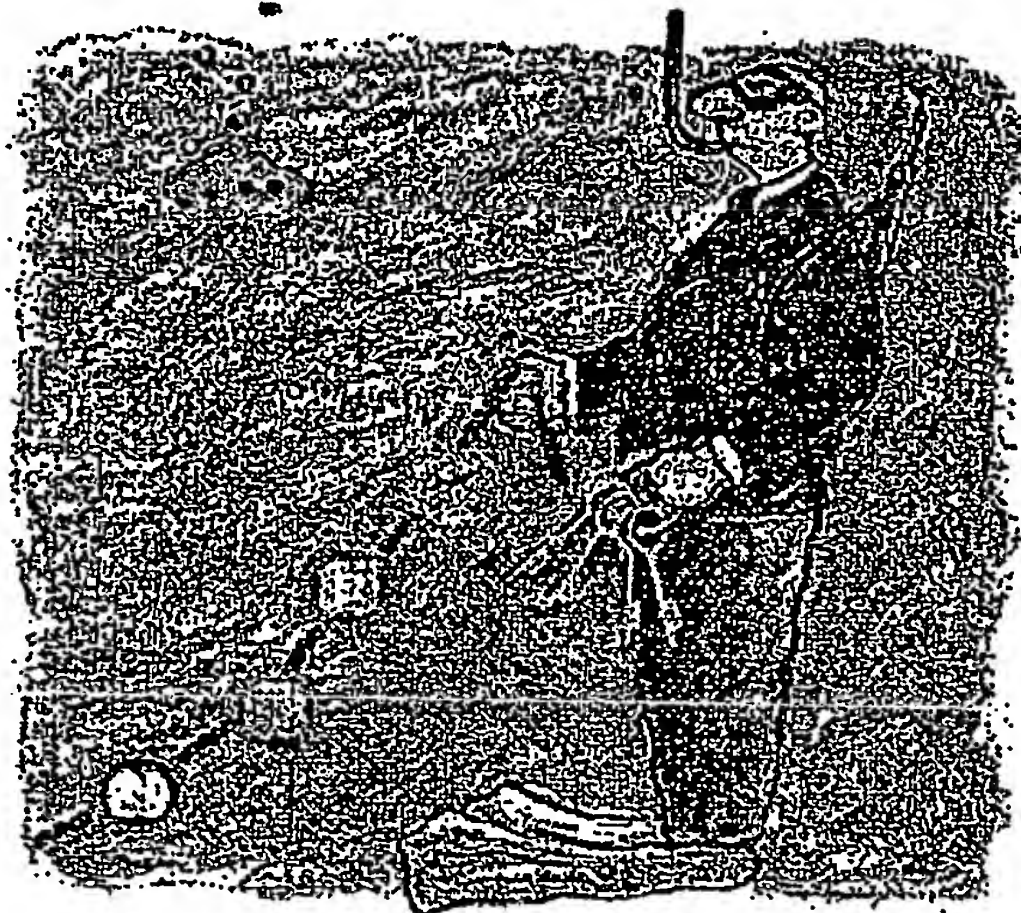


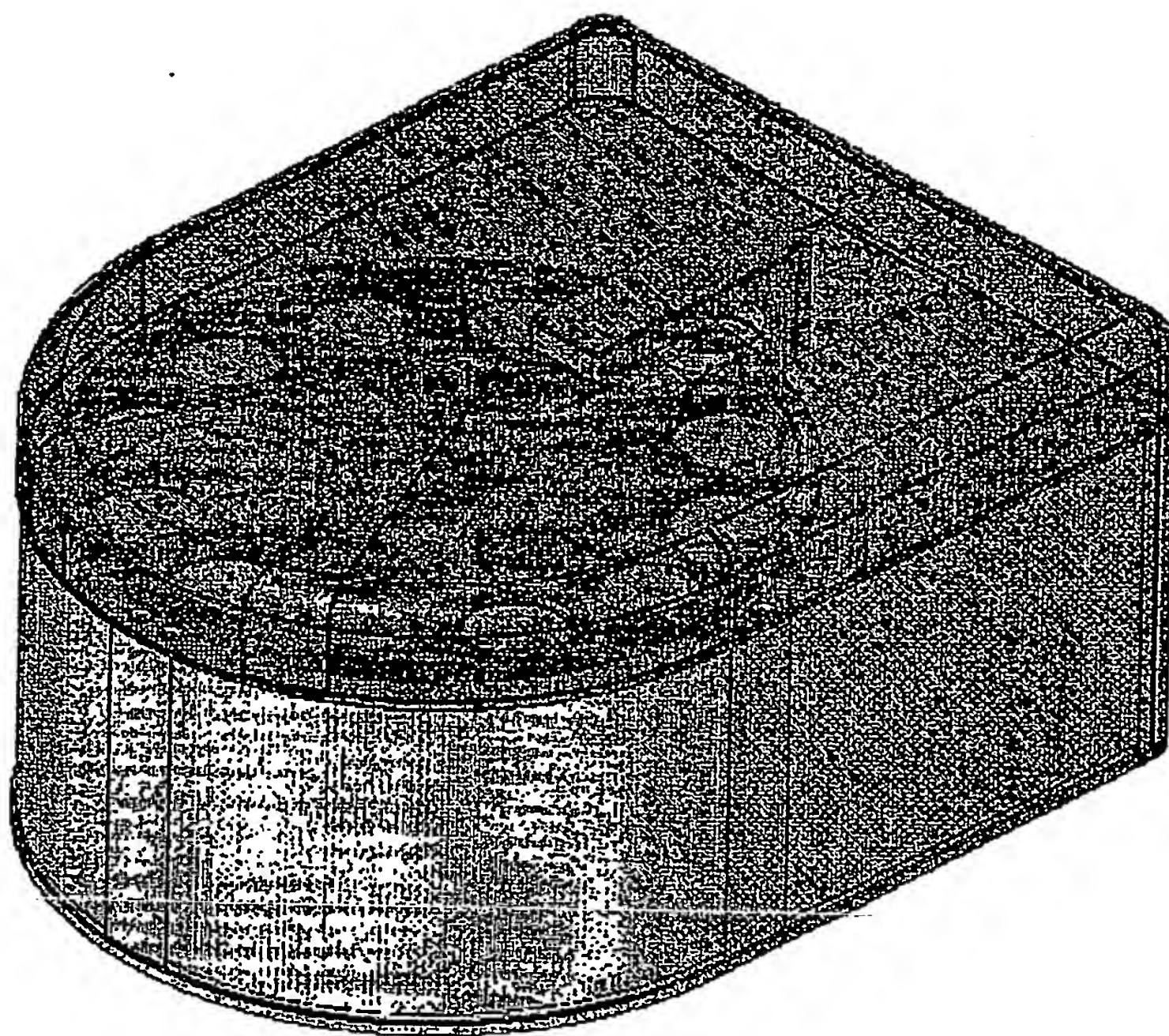

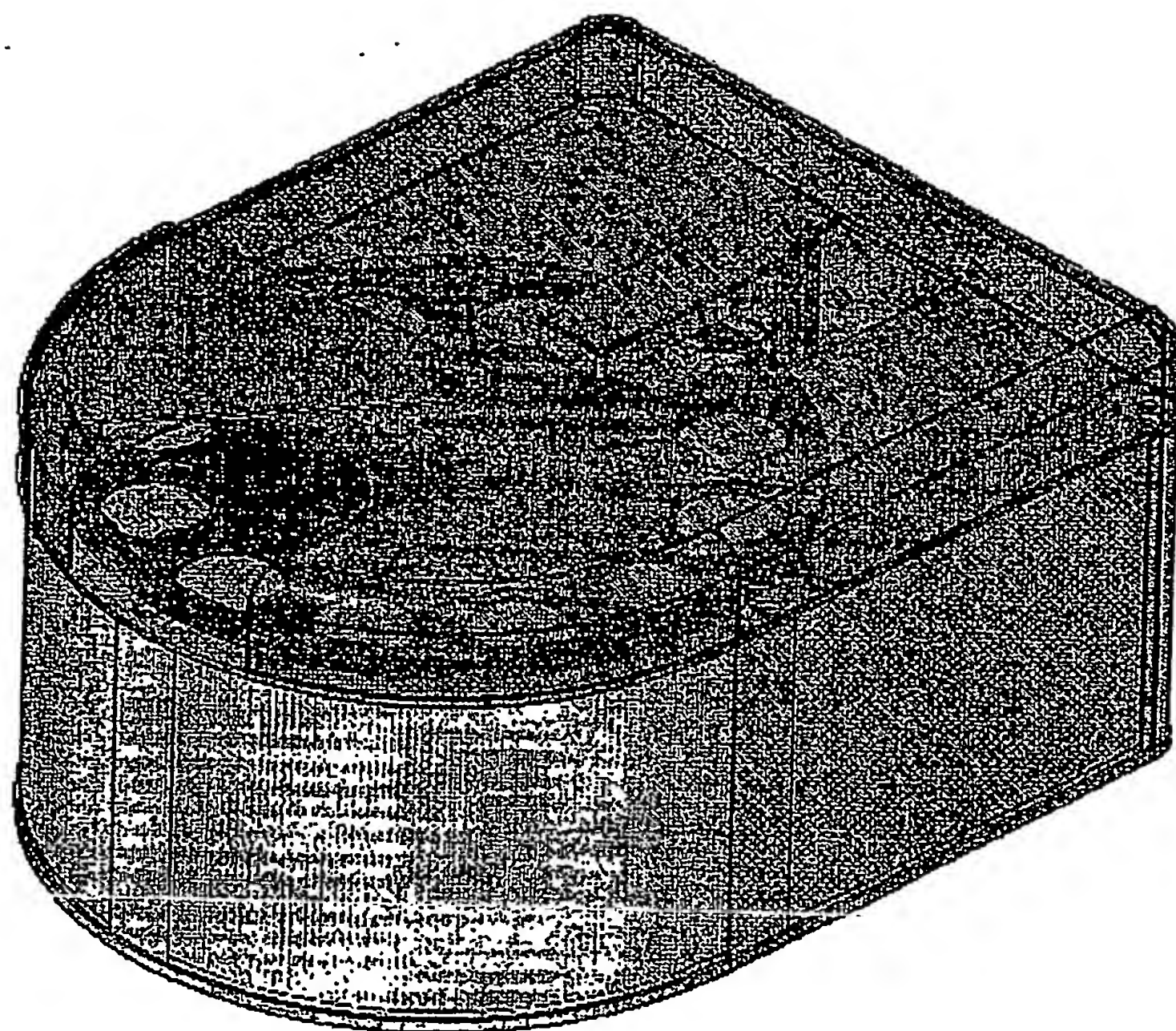





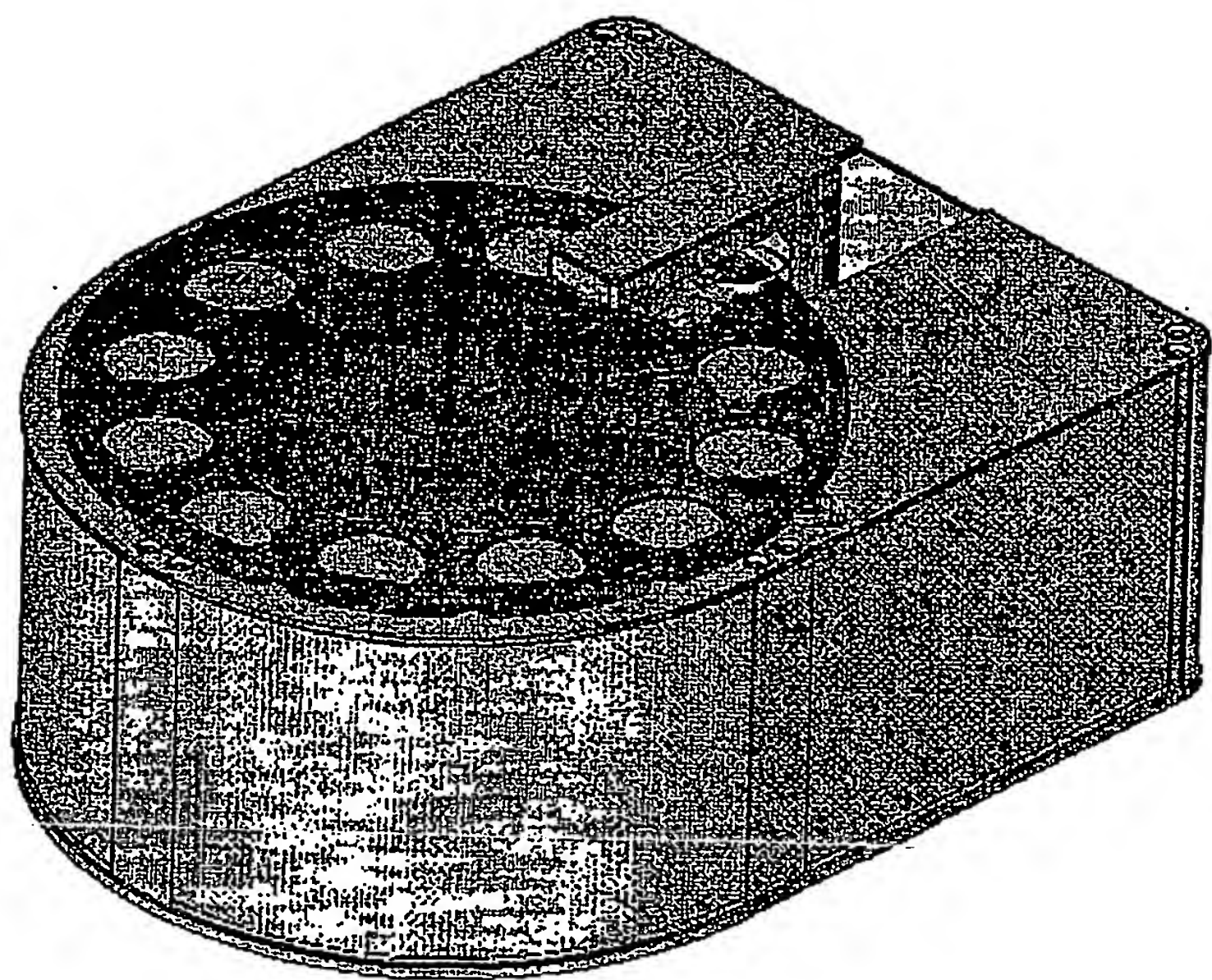

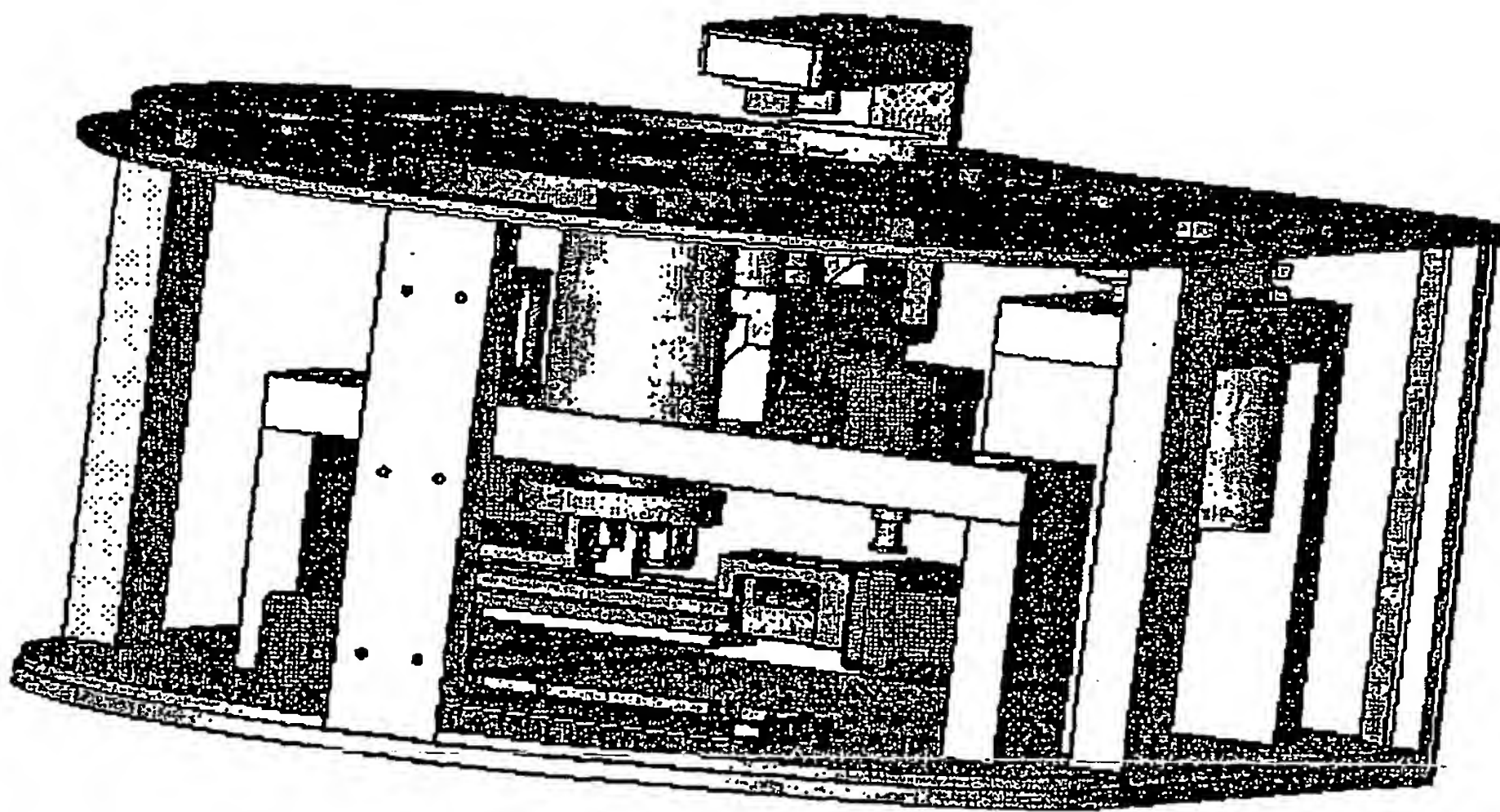

[illegible]

Division of Studies

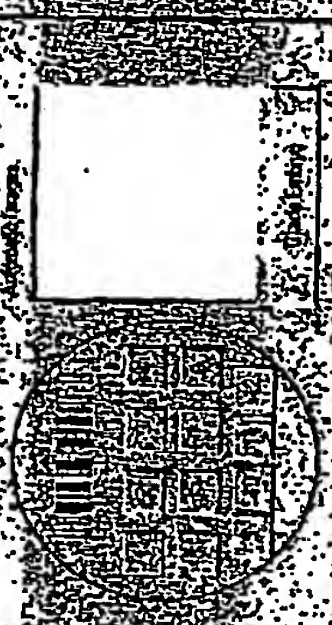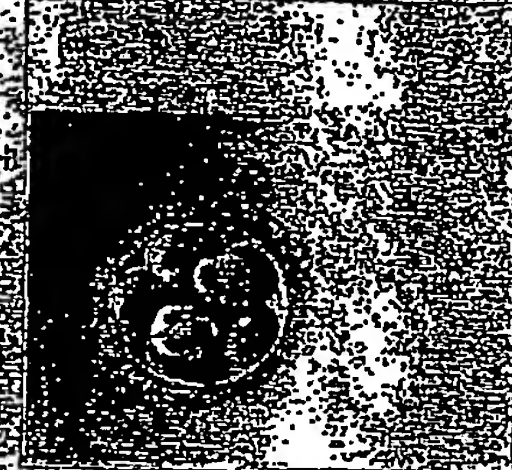

Grade	Registration	Year	Class	Technique	Platform
Nil	-	5/15/02	102-21 PM	Technique 3	1 C-Entrypoint
Cable	-	5/16/02	102-21 AM	Technique 3	1 C-Entrypoint
Cable	-	5/19/02	102-21 PM	Technique 3	1 C-Entrypoint
N	-	5/19/02	102-21 PM	Technique 3	1 C-Entrypoint

[illegible]

The image shows a document page with a grid pattern. The page is heavily obscured by numerous small, dark, rectangular redaction marks. The grid lines are visible, and the redactions are scattered across the entire page, with some larger, more prominent ones near the top and bottom edges.

Unit	Unit 1	Unit 2	Unit 3	Unit 4	Unit 5	Unit 6	Unit 7	Unit 8	Unit 9	Unit 10	Unit 11	Unit 12	Unit 13	Unit 14	Unit 15	Unit 16	Unit 17	Unit 18	Unit 19	Unit 20	Unit 21	Unit 22	Unit 23	Unit 24	Unit 25	Unit 26	Unit 27	Unit 28	Unit 29	Unit 30	Unit 31	Unit 32	Unit 33	Unit 34	Unit 35	Unit 36	Unit 37	Unit 38	Unit 39	Unit 40	Unit 41	Unit 42	Unit 43	Unit 44	Unit 45	Unit 46	Unit 47	Unit 48	Unit 49	Unit 50	Unit 51	Unit 52	Unit 53	Unit 54	Unit 55	Unit 56	Unit 57	Unit 58	Unit 59	Unit 60	Unit 61	Unit 62	Unit 63	Unit 64	Unit 65	Unit 66	Unit 67	Unit 68	Unit 69	Unit 70	Unit 71	Unit 72	Unit 73	Unit 74	Unit 75	Unit 76	Unit 77	Unit 78	Unit 79	Unit 80	Unit 81	Unit 82	Unit 83	Unit 84	Unit 85	Unit 86	Unit 87	Unit 88	Unit 89	Unit 90	Unit 91	Unit 92	Unit 93	Unit 94	Unit 95	Unit 96	Unit 97	Unit 98	Unit 99	Unit 100
Unit 1	Unit 2	Unit 3	Unit 4	Unit 5	Unit 6	Unit 7	Unit 8	Unit 9	Unit 10	Unit 11	Unit 12	Unit 13	Unit 14	Unit 15	Unit 16	Unit 17	Unit 18	Unit 19	Unit 20	Unit 21	Unit 22	Unit 23	Unit 24	Unit 25	Unit 26	Unit 27	Unit 28	Unit 29	Unit 30	Unit 31	Unit 32	Unit 33	Unit 34	Unit 35	Unit 36	Unit 37	Unit 38	Unit 39	Unit 40	Unit 41	Unit 42	Unit 43	Unit 44	Unit 45	Unit 46	Unit 47	Unit 48	Unit 49	Unit 50	Unit 51	Unit 52	Unit 53	Unit 54	Unit 55	Unit 56	Unit 57	Unit 58	Unit 59	Unit 60	Unit 61	Unit 62	Unit 63	Unit 64	Unit 65	Unit 66	Unit 67	Unit 68	Unit 69	Unit 70	Unit 71	Unit 72	Unit 73	Unit 74	Unit 75	Unit 76	Unit 77	Unit 78	Unit 79	Unit 80	Unit 81	Unit 82	Unit 83	Unit 84	Unit 85	Unit 86	Unit 87	Unit 88	Unit 89	Unit 90	Unit 91	Unit 92	Unit 93	Unit 94	Unit 95	Unit 96	Unit 97	Unit 98	Unit 99	Unit 100	



[illegible]

[illegible]

Is Get new Dist

Technician 2	Dist 5h	12456789	Remarks
12456789			
Get new Dist			



Drop Property

6/16/02 2:32:23 PM Technician 3

6/17/02 8:32:23 AM Technician 3

6/18/02 1:32:23 AM Technician 3

6/19/02 12:32:23 AM Technician 3

6/18/02 11:32:23 PM Technician 2

2 Mil

2 ZPH

9:3 Cells

2:7 Cells

2:8 Cells

Start Time

6/18/02 2:32:23 PM

6/19/02 11:32:23 PM

End Time

6/18/02 11:32:23 PM

6/19/02 11:32:23 PM

Tested

3

2

9 Cells

Technician 2

Technician 3

[illegible]



#### *5.1.1.1. Bovine Semen*

The worldwide use of bovine semen in animal husbandry is enormous, with approximately 100 million doses administered each year. With costs per dose ranging between approximately \$4-50/dose, and with a mean of about \$10/dose, this total industry has a value of about \$ 1 billion per annum. IMT estimates that at a charge of \$0.90 per dose, it's market share would be about \$100 million per annum.

Use of IMT's MTG™ 525 equipment in the freezing of bull semen has been shown to yield a greater proportion of viable sperm than other conventional methods. This leads to more doses of semen available per ejaculate, a major commercial advantage for the user. In addition, because the IMT technology reduces significantly the damage to sperm cells, it may be possible to achieve successful insemination with smaller doses of semen, an additional advantage of this technology. Field trials are now underway at Cogent Breeding Ltd (UK) to demonstrate this added advantage. Additional potential advantages of IMT technology – the ability to freeze larger volumes of semen in a single container, the ability to refreeze semen specimens after initial thawing, and the freezing of semen specimens that have been sex – sorted – are some additional applications that are being investigated by IMT. All these variations are desirable additional applications of the IMT technology that would have significant operational and commercial benefit to the user. By developing and demonstrating these applications, IMT will be able to secure its place in the market for a long time.

IMT currently has the MTG™ 525 module in routine use in the UK (Cogent). It is also under evaluation in Switzerland. IMT anticipates that an additional 5-6 machines will be placed during the second half of the year 2002.

#### *5.1.2 Equine Semen*

Equine semen freezing is not as far advanced in the world market as is bovine semen, but the worldwide market is growing rapidly. Stallion semen is more difficult to freeze than bovine semen, and the semen of about 25% of stallions seem to be resistant to freezing. Success rates for equine artificial insemination is about 40-50%.

As currently practiced, stallions selected for siring are often shipped at sometimes great distances so that natural insemination can take place. This approach is costly and inefficient. In contrast, artificial insemination permits the "banking" of selected stallion semen. Insemination of a mare can take place at the convenience of the owners, even long after the stallion has died. Thus, the use of semen storage offers the horse's owner a form of insurance over his (usually large) investment.

Preliminary studies indicate that IMT freezing technology can improve semen recovery by about 20%. More importantly, semen from stallions that could not be frozen previously has been successfully frozen by this technology and has resulted in pregnancies. This opens up a large and remunerative field for IMT's technology. Based on a \$30/dose charge, IMT estimates this current annual market at about \$300 million. IMT, in collaboration with Cogent is now continuing to perform field trials in Europe to further document the efficacy of its equipment in the freezing of equine semen. Results are anticipated by the end of year 2002.

#### **5.1.3 Porcine Semen**

IMT has recently added the freezing of porcine semen to its growing list of animal artificial insemination applications. Studies in IMT's laboratories indicate that, when using the MTG™ 525 module, there is an excellent, 95% recovery of viable sperm. Field trials in which porcine semen frozen with IMT's technology will be used to inseminate and produce viable offspring are scheduled to commence in the Summer of 2002.

Average production cost of one insemination dose is about \$3.5 when produce in the farm (disregarding the boar's genetics). An average farm produces about 15,000 doses annually. This is a market that is in its relatively early stages of development, and its size is difficult to estimate. IMT anticipates that an MTG™ technology will change this field and that there will be newly established artificial insemination farms supplying frozen porcine semen. IMT estimates that it will be able to charge a royalty fee of \$0.90 per dose of porcine semen frozen.

## Successful pregnancies in cows following double freezing of a large volume of semen

### Abstract

The objective of the following paper is to describe a new technology for large volume and double freezing of semen in 12 ml test tube. Semen from two different bulls was frozen with a new technique using 12 ml test tube and was refrozen after thawing in mini straws. All freezing was done in a "Multi thermal gradient" (MTG) freezing apparatus, which moves the container at a constant velocity (V) through a thermal gradient (G) producing a controlled cooling rate  $B = (G) \times (V)$ . Each of the two bulls ejaculated were evaluated for post thaw motility in the lab and then in a field trial which was carried out in a split sample mode. We inseminated 105 cows after double freezing/thawing cycle, and another 123 cows were inseminated with semen frozen in mini-straws and a conventional method. Results showed a  $75 \pm 5\%$  post thaw motility after freezing a 12ml test tube and  $50 \pm 5\%$  after second freezing/thawing in mini-straws, respectively. Controlled vapour freezing showed a  $60 \pm 10\%$  post thaw motility. Results of the field trial showed a pregnancy rate of 44% (47/105) for the double freezing group in comparison to 45.5% (56/123) for the controlled group. These results can be beneficial for large volume freezing, and therefore for bull semen cryobanking in a large volume which will be followed by second freezing in a regular insemination volume.

### Bull semen Cryobanking

Cryobanking of semen has had a major impact on dairy cattle genetic breeding. In addition to its role in young bull genetic breeding, cryobanking of bull semen is an important backup for sufficient insemination doses in cases of disease, infertility or mortality.

Freezing and storage of semen is done regularly using mini ( $\frac{1}{4}$ cc) or midi ( $\frac{1}{2}$  cc) straws. However, cryobanking of a large number of straws is time consuming, expensive and requires a lot of storage space and liquid nitrogen. An alternative procedure which will reduce these expenses could be the freezing of a whole ejaculate in one test tube (12ml) and only when needed (when the bull is a "proven bull") the test tube will be thawed and then be refrozen in regular mini straws. We describe here the use of a new technology for large volume (whole ejaculate) freezing/thawing and refreezing in mini-straws.

### MTG technology

Our novel freezing technology is based on "Multi-thermal gradient (MTG<sup>®</sup>, IMT, Israel) (1)" directional solidification and is used mainly for freezing sperm and large



tissue. The semen in the test tube is moved at a constant velocity ( $V$ ) through a linear temperature gradient ( $G$ ) so the cooling rate ( $G \times V$ ) and ice front propagation are precisely controlled (Fig. 1).

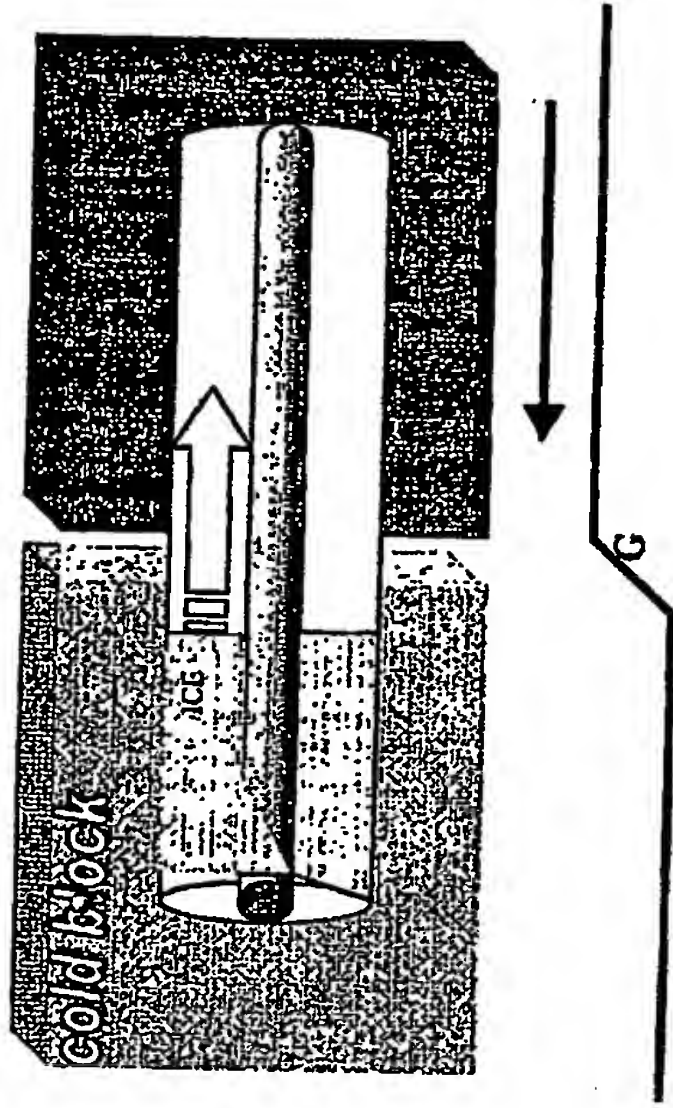


Fig. 1 Schematic design of the MTG freezing

This method also enables the incorporation of controlled seeding into the freezing process. When any liquid is cooled below its freezing point, it remains a liquid, in an unstable super-cooled state, until freezing starts at randomly distributed nucleation sites and spreads throughout the entire volume of the liquid. As discussed above, in the conventional equiaxed method of freezing, ice grows with uncontrolled velocity and morphology, and may disrupt and kill the cells of the samples. Ideally, the velocity of the freezing front should be such that the ice morphology does not disrupt the cells or tissue. However, the rate of cooling appropriate for favorable ice morphology may not be appropriate for other desired outcomes of a sample's freezing protocol. The laterally varying gradient used in our technology allows cooling to proceed at differing rates under varied temperature regimes, thereby facilitating full control over nucleation and ice crystal morphology. This technique allows very precise control of the cooling rate (0.01 to 1000°C/minute) within a large volume.

The freezing apparatus can control ice crystal propagation by changing the thermal gradient ( $G$ ) or the liquid-ice interface velocity ( $V$ ) and so optimizing the ice crystal morphology during freezing of cells and tissue. The rate of cooling also affects the

morphology of the intercellular ice crystals (3): morphologies such as closely packed needles kill cells by external mechanical damage (unpublished observation). Thus, maximizing the survival rate of cells subjected to freezing and thawing requires careful control of the freezing process i.e. interface velocity. Using a cryomicroscopy observation we found that survival of sperm shows biphasic curve where at a very slow velocity ice will grow in a planar form which will kill all cells. At higher velocity ice crystals will form secondary branches and survival will increase, however at higher velocity (i.e.  $300\mu\text{m/sec}$ ) ice will start to form "needle-like" ice crystals which will decrease PTM, but in a higher velocity will permit very high survival (fig. 2) depending on the space between the ice crystals (4). Finally, at very high velocity (i.e.  $>3000\mu\text{m/sec}$ ), directional solidification will not occur and survival will decrease.

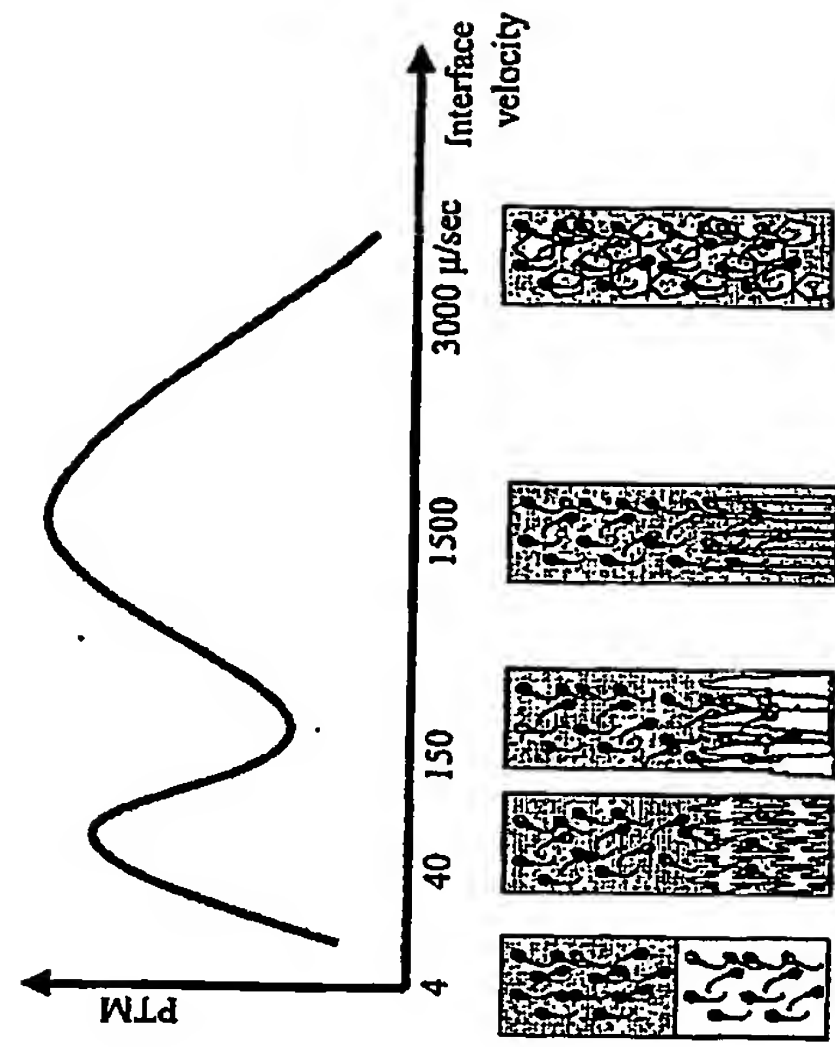


Fig. 2 Effect of interface velocity (V) on ice crystals morphology and sperm post thaw motility (PTM)

Heat transfer problems associated with large volume freezing

In a conventional slow-freezing method, temperature of the chamber is dropped in a controlled stepwise manner. This method is based on using multidirectional (equiaxed) heat transfer to achieve a rate of temperature change in the sample that depends on the thermal conductivity and geometrical shape of the container and of the biological material within it (5). The thermal gradient within the sample is determined implicitly by the temperature of the chamber and the thermal conductivity of the materials of the sample, and is not directly controllable. Furthermore, the ambient temperature gradients within the freezing chamber and the unreliability of temperature recording measurements (6) add to the difficulty of achieving the optimal cooling rate in a large volume sample.

#### Cryobanking of large volume semen

Each of the ejaculate was tested for semen concentration and motility (>70%) before dilution. We used AndroMed<sup>®</sup> (minitub, Hauptstrabe, Germany) for the semen dilution to have a final concentration of  $1.5 \times 10^6$  sperm/ml.

Freezing of a whole ejaculate was done in a special test tube (12ml) in which the central part is a hollow channel. Heat transfer is opposite to the test tube movement and is parallel to the tube length axis (fig. 1). The empty channel in the middle of the large test tube facilitate directional freezing and rapid thawing in the inner side of the test tube.

Sperm PTM after freezing in a large volume was very high. We found a survival rate of 90-100% (normalized PTM) in the two bulls we cryopreserved in the MTG technique. These results were superior to MTG freezing using mini straws (data not shown), which suggest the benefit of using MTG freezing of large volume for sperm cryopreservation. Results shows a  $75 \pm 5\%$  post thaw motility after freezing a 12ml test tube and  $50 \pm 5\%$  after second freezing/thawing in a mini-straws, respectively. Controlled vapour freezing showed a  $60 \pm 10\%$  post thaw motility which were lower then the results after MTG freezing of mini straws.

The large volume freezing may be very useful for cryobanking of bull semen, for example, AI centre that have a bank of 10,000 straws which are made from 25 ejaculates (400 straws /ejaculate). We calculated that these 10,000 straws will fit into 13 goblets (750 straw/goblet). In comparison, when we freeze a large volume (12ml test tube) the 25 ejaculates will be frozen in 25 test tubes which will be stored only in 2 goblets. This means that we need 6.5 time more goblets using straws in comparison to test tube freezing. In this case, the present method gives a capability to have a bank of "waiting bulls" in some of the AI centers which presently do not use a semen cryobanking. In addition this method will save money in labour and consumables (filling, printing, LN for freezing and for storage etc.).

In conclusion, the MTG technique could be very useful for large volume cryopreservation and double freezing for sperm cryobanking.

#### References

1. Arav A., Device and methods for multigradient directional cooling and warming of biological samples. US Patent (1999) 5,873,254

2. Foote R., and Park J. Factors affecting preservation and fertility of bull sperm: a brief review. *Reprod. Fertil. Dev.* 5, (1993) 665-673
3. Karlsson J. O. M., and M. Toner., Long term storage of tissues by cryopreservation: critical issues *Biomaterials* 17 (1996) 243-256
4. Koebe H. G., Werner A., Lange V., Schildberg F. W., Temperature gradients in freezing chambers of rate-controlled cooling machines. *Cryobiology* 30 (1993) 349-353
5. Mazur P., The role of intracellular freezing in the death of cells cooled at supraoptimal rates. *Cryobiology* 14: (1977) 251-272
6. Mazur P., Cole K.W., Role of unfrozen fraction, salt concentration, and changes in cell volume in the survival of frozen human erythrocytes *Cryobiology* (1)(1989) 1-29



multigradient notes

## Multigradient

Notes from interview with Amir Arav, 26 June 27-06-2002

Based on US 5,873,254.

Also based on the cryopreservation part of US 60/345,643.

Equine: most stallion sperm can't be frozen. We get higher recovery than others.

Largest volume of semen frozen in prior art in test tube = 5 ml. (Larger volume of pig semen has been frozen in bags.)

Our innovation: large volume freezing.

Horse prior art: 0.5 ml. One insemination needs four to eight 0.5 ml semen samples.

### Why large volume?

For example, for breeding cattle. Genetic breeding. Test young bulls for production. Collect semen. Inseminate heifers. Heifers give birth to heifers. Inseminate the daughters. See how much milk the daughters give. This takes 4 to 5 years. Only one out of every 14 candidate bulls is selected and something may happen to the best bull during the 4-5 years, like the bull may die. Therefore, need to put semen in bank. Store 10,000 to 50,000 0.25 ml straws per bull. It takes 25 days to collect 10,000 straws, and lots of liquid nitrogen for storage.

For volume of one ejaculate: We do one test tube, 12 ml sample, 1:2 dilution instead of 400 to 600 straws, 1:10 dilution. At the end of the 4-5 years, the selected bull's semen is thawed and refrozen in regular straws.

Prior art concept: can't freeze samples bigger than 0.5 ml, can't freeze concentrated sperm.

Prior art: 50,000,000 sperm cells per ml.

Us: 500,000,000 sperm cells per ml.

Equine and boar semen freezing without centrifugation.

Prior art: do centrifugation of semen to remove seminal plasma before freezing.

Sperm concentration in semen is low in stallion and boar: 50,000,000 to 600,000,000 sperm cells per ml.

To remove plasma: centrifuge and wash the sperm. This damages the sperm. Then add extender to get concentration needed for insemination: 1,000,000,000 to 6,000,000,000 sperm cells total.

Us: dilute the plasma. Go down to 20,000,000 sperm cells per ml. Then need 50 ml for one insemination. Large volume freezing allows freezing one insemination (10 ml x 5 or 50 ml x 1) at once.

How?

Multigradient freezing of rotating test tube. Rotate the tube around its longitudinal axis during the freezing.

multigradient notes

Advantage: mix the solution in front of the ice front. This dilutes the concentration of salt being expelled from the ice.

Rotating tube also used for warming.

Rotating tube also used for freezing partly filled test tube. Spread solution in annulus to get high surface area.

Rotating the tube gives better thermal contact between the solution and the metal heat exchanger through the walls of the test tube.

The thermal contact with the block is always best on the bottom side of the tube, and that thermal contact gets spread around when the tube is rotated.

If the tube is partly filled with a sample, you wind up with frozen sample on the wall of the tube and air along the axis of the tube. That's better if then you want to lyophilize (freeze dry) the sample.

Rotating the tube keeps the sample mixed and homogeneous during freezing.

Alternative: hollow (double walled) test tube.

Piston sits inside central channel of test tube and removes heat from the central channel. More efficient heat exchange during cooling.

For heating: put hollow test tube in water bath so circulating water flows through the central channel. Warming rate in center of tube and on outside of tube is the same.

Unlike US 5,873,254 seeding, need to inject liquid nitrogen at bottom of test tube. A section at the base of the test tube is arranged to exclude sperm but include liquid, for example by putting glass balls as liquid trap in base of test tube, so only the liquid is frozen for seeding.

Hollow tube can be glass or plastic. Hollow tube need not rotate.

Hollow tube also has a roughened section on side for manual marking.

Warming or thawing small straws (0.25 cc or 0.5 cc samples).

Faster is better, but can't go fast from liquid nitrogen temperatures or heat stress will crack the sample.

OTOH, overheating leads to denaturation, and warm cryoprotectant can damage sperm.

Faster is better to prevent recrystallization at  $-10^{\circ}\text{C}$ .

Machine has one block at uniform high temperature:  $38^{\circ}\text{C}$  to  $100^{\circ}\text{C}$ , with  $90^{\circ}\text{C}$  being preferred. Put straw in hole through block. Move straw through the block via the hole at constant velocity. 6mm/sec optimal. Block is 2 cm thick so 3.3333 second duration. Then out to ambient air.

From when straw leaves liquid nitrogen to when straw emerges from block should be less than 50 seconds, preferably about 30 seconds.

Pull straw from liquid nitrogen: goes up to  $-30^{\circ}\text{C}$  just by being in ambient air. Put straw through block: in 3 seconds, go to room temperature.

Another way to warm the test tube: like the straws.

Another way to warm a prior art test tube:

multigradient notes

Plunge into a water bath warmer than 37°C, with 70°C preferred. So no contact between ice and walls. Then drop contents of test tube into high volume pre-warmed dilution (insemination) solution. Mix.

Alternative: kept the test tube in hot water while "stirring" with the test tube to get uniform thawing.

Freezing and thawing test tube with controlled rolling (rotating) system

Inventors: Arav Amir, Meir Uri

We developed a device which has controlled rolling system of round container (i.e. test tube) during freezing and thawing.


The advantageous of this device are:

1. A better heat transfer between the container and the copper blocks.
2. A formation of air bubble in the centre of the container/ a thin layer on the wall of the container.
3. A controlled propagation of ice crystals in parallel to the container wall.
4. A continuous mixing of the solution during the freezing and thawing.
5. A preparing of large surface for the purpose of sublimation for freeze drying.





CROSS - SECTION OF HOLLOW TEST TUBE

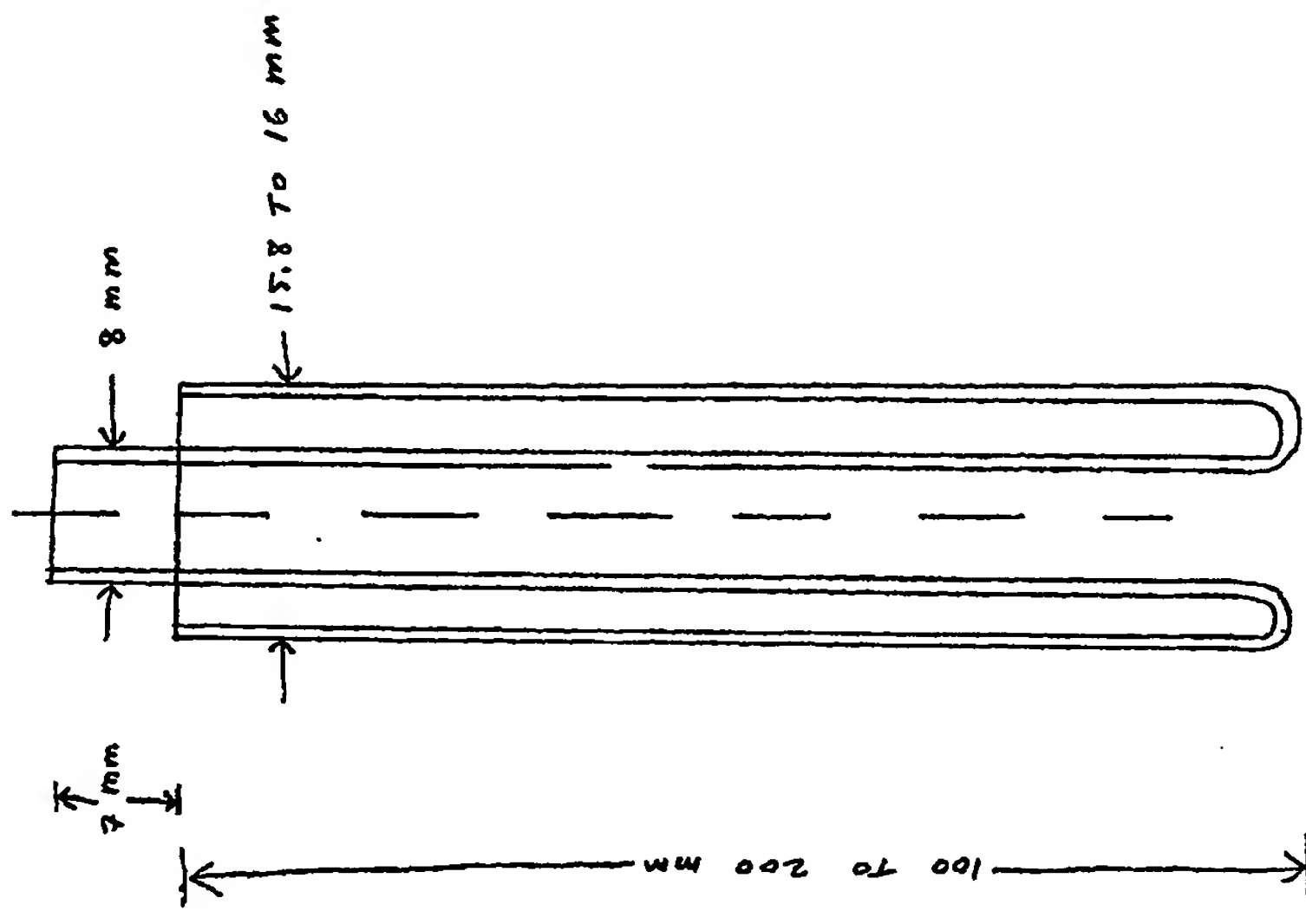

COGENT BEACHIN STLD

DOUGLAS FRIEDMAN FIRST EXPERIMENT

The first experiment was to find the best dilution rate for the bacterium (1:1, 1:2, 1:3) for 12 ml. Tubes for dilution 1:1 the serum was diluted with an extender of 10% glycerol for a final concentration of 8% glycerol. For dilution 1:2 the serum was diluted with an extender of 17.5% glycerol for a final concentration of 8% glycerol. For dilution 1:3 the serum was diluted with an extender of 25% glycerol for a final concentration of 8% glycerol.

FRIEDMAN

The serum was taken from two bulls and divided in three parts for the different dilution rates, the same was frozen in the NITRO-310 machine of 12 ml. at a velocity of 1.0 mm/sec. The temperature for block A was 60 and for block B -80C. normal serum was frozen in Opti Cool machine as control.

THAWING

The tubes were held in the air for 140 sec. and then submerged in a water bath at 40C for 25-30 sec. the contents of the tube was added to a freezing solution at 37, the serum quality was assessed in the microscop and then cooled down to 20C, the second freezing was done in 0.5 ml. serum for AUTO-550 and 0.25 ml. serum for Opti Cool machine as control.

RESULTS

First freezing results

Case	bull	1:1 ratio	1:2 ratio	1:3 ratio	normal serum
02.04.02	Almas	49.9	63.3	63.3	Control D.C.
	C. hakeam	60.0	70.0	47.5	70.0
AVG		55.3	66.7	55.4	60.0

Second freezing results

	1:1 ratio	1:2 ratio	1:3 ratio	Control D.C.
Almas	20.0	48.8	25.0	25.0
C. hakeam	30.0	33.0	30.0	30.0
AVG	25.0	40.9	27.5	27.5

Conclusion: The results show the best dilution rate for freezing in 12 ml. tubes is 1:2 Serum-Accelerator 60.7% and lower results were found for 1:1 and 1:3, 63.3% and 60.4% respectively, the results for the C. Cool machine (control) was 65%. In the second freezing the best result was obtained in 1:2 ratio

DOUBLE FREEZING SECOND EXPERIMENT.

After the first experiment we find that the best results for freezing 12 ml. Tube was 12 semen/ender dilution. The second step was to check a big number of bulls based on the results of the first experiment. for better comparison in this experiment we used 0.215 ml. straws with the same sperm concentration 20X10<sup>6</sup> for all treatments.

MORTALITY %

	C. Ryan	L. Wilam	H. Concoro	Alasia	L. Lucan	W. Glunnen	G. Dalton	AVG
Raw Semen	70	80+	80	90	90+	80	80	78.0
After first freezing tubes	60	65	60	65	60	60	60	61.4
After first freezing D. Cool (control)	65	65	60	70	60	70	60	70.0
After second freezing straws IMT	28	50	18	50	60	40	50	40.9
After second freezing straws D. cool	15	45	40	40	40	60	40	34.5

Conclusions: the average of the results show

- 1) After first freezing Digi Cool machine it better than first freezing in 12 ml. Tubes in IMT machine 70% vs 61.4%.
- 2) After second freezing IMT machine shows a little better results the Digi Cool machine 40.9% vs 34.5%
- 3) There are individual differences between the bulls in the freezing capacity of the semen and between the two freezing techniques, for example the bull Concoro freeze better in D. cool machine than IMT machine 40% vs 18%
- 4) We can conclude that is possible to perform a field trial after second freezing with both freezing techniques in the second freezing, with those bulls that it semen quality varies from 40% to 80% in D. Cool and IMT machine in the second freezing the best bulls for inseminations trials are: L. Wilam, Alasia, L. Lucan, W. Glunnen and G. Dalton

SCORED SEMEN MOTILITY, DEADLINE RATIO AND MEMBRANE INTEGRITY AFTER THAWING

The test experiment:  
The aim of this experiment was to test the best thawing time for 0.5 ml straws at 90 degrees Celsius.  
For straws of 0.25 ml was found in previous work to be 2 seconds in 90 degrees Celsius the best thawing time.  
At the beginning of each experiment were tested two, three or four straws of 0.5 ml and thawed at different thawing times. Two test were applied for the assessment of the semen:  
1) Visual motility after 5 min. of incubation at 37C  
2) Propidium iodide fluorescent test for dead/live ratio dead cells red colour and live cells blue colour.

The results:

Motility of semen at different thawing times of 0.5 ml straws at 90C

date	bull	3	4	5	6
24.04.02	Meadowbrook			60	30
30.04.02	Shaker		60	65	
31.04.02	Lucky	60	60	60	
01.05.02	Shaker	35	45	60	
02.05.02	Booker	80	40	75	60
05.05.02	Principal	20	30	60	60
07.05.02	Lucky		60	70	60
10.05.02	Courier		60	70	60
17.05.02	Shaker		60	60	40
18.05.02	Courier		30	60	60
AVG		41.3	47.2	61.0	47.1

Propidium iodide test dead/live sperm ratio

date	bull	3	4	5	6
24.04.02	Meadowbrook			75	67
30.04.02	Shaker		61	65	
31.04.02	Lucky	55	64	65	
01.05.02	Shaker	63	73	73	
02.05.02	Booker	61	79	80	71
05.05.02	Principal	51	68	77	78
07.05.02	Lucky		68	61	60
10.05.02	Courier		68	64	66
17.05.02	Shaker		77	63	62
18.05.02	Courier		63	72	63
AVG		62.6	68.3	70.5	69.3

Results: The average of the results show a better motility and Propidium iodide test 68% and 70.5% respectively after five sec. thawing in water temperature at 90C for 0.5 ml straws.  
Lower results were found in 3, 4 and 6 sec. thawing.  
Conclusion: The best results for both test were five sec. thawing at 90C. for 0.5 ml straws.



**SHORTED BEMEN MOTILITY, DEADLYIVE RATIO AND MEMBRANE INTEGRITY AFTER THAWING**

The second experiment:

Based on the first experiment the content of bechre bulb that normally produces semen for cooling was tested in two temperatures 37C and 90C and two volumes of straw 0.25 ml and 0.5 ml in order to find the best method for freezing and thawing semen in the MTG-623 and MTG-550 machines.

The freezing time of 1min at temperature of 37C was the same for 0.25 ml straw and the thawing time for temperature of 90C was 2 sec for 0.25 ml straw and 5 sec for 0.5 ml straw.

The freezing velocity for 0.5 ml straw was 1.5 mm/sec. With a seeding time of 60 sec and for the 0.25 ml straw the freezing velocity was 2.0 mm/sec. With a seeding time of 11 sec.

three test were applied for the assessment of the semen quality.

1) Visual motility after 5 min. incubation

2) Propidium iodide test for dead/live ratio

3) OAT test for membrane integrity, the sperm is incubated in a hypotonic solution for 40 min with a fluorescent dye that can identify the attached sperm, we used propidium iodide stain. then is stained

The results:

date	bulb	Motility after thawing			
		0.25 ml 37C.1min	0.25 ml 90C 2sec	0.50 ml 37C.1min	0.50 ml 90C. 5sec
23.05.02	Dazzler	60	60	60	60
24.05.02	Meadowstream	50	60	50	55
26.05.02	N. Omar	65	60	50	60
30.05.02	Shakar	60	60	60	60
31.05.02	Lucky	65	60	70	60
01.06.02	Shakar	60	55	60	60
02.06.02	Booker	60	60	60	60
05.06.02	Principal	65	60	60	70
07.06.02	Lucky	55	65	63	73
10.06.02	Courier	60	70	50	70
17.06.02	Shakar	60	70	70	65
18.06.02	Courier	60	60	60	65
AVG		57.8	62.0	57.8	65.3

**Propidium iodide dead/live ratio**

date	bulb	Propidium iodide dead/live ratio			
		0.25 ml 37C.1min	0.25 ml 90C 2sec	0.50 ml 37C.1min	0.50 ml 90C. 6sec
23.05.02	Dazzler	69	68	67	70
24.05.02	Meadowstream	62	66	64	74
26.05.02	N. Omar	43	49	75	63
30.05.02	Shakar	54	69	69	69
31.05.02	Lucky	67	60	63	55
01.06.02	Shakar	64	66	78	74
02.06.02	Booker	73	70	67	79
05.06.02	Principal	64	62	66	62
07.06.02	Lucky	79	77	33	70
10.06.02	Courier	77	76	68	76
17.06.02	Shakar	61	69	69	68
18.06.02	Courier	62.3	64.8	60.7	69.3
AVG					

## CORT 40 min. incubation (membrane integrity test)

	0.25 ml. 37C. 1min	0.25 ml. 90C. 2min	0.50 ml. 37C. 1min	0.50 ml. 90C. 5min
23.04.02 Dazler	43		55	
24.04.02 Moudystream	31	37	34	43
26.04.02 N. Omar	29	35	45	63
30.04.02 Shahr	40	42	61	57
31.03.02 Lucky	57	44	32	35
01.05.02 Shahr	58	67	44	38
02.04.02 Booker	43	61	35	65
05.05.02 Principal	50	38	58	46
07.04.02 Lucky	35	63	30	43
10.04.02 Courier	55	54	34	51
17.05.02 Shahr	70	63	62	70
19.04.02 Courier	28	44	28	41
AVG	43.8	48.0	39.8	48.3

The results: The best results in molality were for the high thawing temperatures (90C) for the volumes 0.25 ml and 0.5 ml, 65% and 65.3% and the lower results were for the low temperatures (37C) for 0.25 ml and 0.5 ml with 57.5% and 67.8% respectively.

For the Propidium iodide test (dead/live ratio) the best results were for 0.25 ml and 0.5 ml strains thawed at high temperatures 90C, 64.8% and 65.3% respectively, and lower results were found for 0.25 ml and 0.5 ml strains thawed at low temperatures 37C, 52.3 and 60.7 respectively

For the CORT test (membrane integrity of the sperm cells) the best results were for 0.25 and 0.5 ml strains thawed at high temperature 90C, 48.0% and 49.3% respectively and lower results were found for 0.25 ml and 0.5 ml thawed at low temperatures 37C, 48.0% and 39.8% respectively.

## Conclusions:

The strains thawed at high temperatures shows an improvement in motility, Propidium iodide, and CORT test means that less damage is caused to the sperm cells when the thawing is doing at high temperatures 90C (fast thawing). The volume of 0.5 ml strains show a slightly better results in all three test than for 0.25 ml strains when thawed at high temperatures 90C, this results show an advantage for big volumes, that seems to give better protection to the sperm cells under drastic changes in temperature, no differences were found for 0.25 ml and 0.5 ml strains when were thawed at low temperatures 37C except for the CORT test, the better results were for 0.25 ml strains than 0.5 ml strains 45.8% vs 39.8% respectively.

7th 7135 (2011) 6

CURRENT RESEARCH STUD

Final tuning, comparison between 0.5 ml. Straws frozen in MTC-550 vs. 0.25 ml. straw frozen in MTC-625 and Digt Cool machine as controls.

The aim of this experiment was to compare 0.5 ml. Straws frozen in MTC-550 in three different velocities 1200, 1600, and 1800 mm/sec. vs. 0.25 ml. at 2000 mm/sec. the best results for this machine. The temperatures for both machines were 5C for the first block and -50C for the second block, the thawing temperatures of the straws was 37C.

The data is from eleven bulls, that give semen for normal freezing. The concentration of the semen was 20X10<sup>6</sup> per straw, the semen was diluted at room temperature and cooling down at the temperature of 5C, the straws were filled and then frozen in the different machines.

Results:

The motility of semen, frozen in MTC-550 at different velocities compared to MTC-625 and Digt Cool machines

Bull	Raw semen motility	V-1200 MTC-550	V-1600 MTC-550	V-1800 MTC-550	V-2000 Digt Cool MTC-625 Machine
Notulon	80	30	50	5	50
Draumer	80	50	75	5	55
C. Laurel	80	40	85	5	60
Palton	80	50	75	2	75
Ladiner	70	50	60	5	60
Pierpo	80	65	70	65	70
Cashlan	70	60	60	6	65
C. Gracklor	80	30	45	5	50
Pushie	40	30	30	2	30
C. Herriot	60	60	50	2	75
S. Nilton	60	50	50	10	40
AVG	72.7	45.0	58.4	9.2	59.9

Conclusions:


The best results were found for the MTC-625 at a velocity of 2000 mm/sec. 60% and for MTC-550 at a velocity of 1500 mm/sec. with 58.4% motility. Lower results were found for Digt Cool machine conventional freezing with 49.5% and for MTC-625 at a velocity of 1200 mm/sec. with 45% motility, the worst results were obtained for MTC-550 at a velocity of 1800 mm/sec. with 9.2% motility.

Bull sperm drawing in drawing machine

#	T-10	Vel/time	Air time	% motil	SVBR	PI	SVBR	PI	SVBR	PI
1	80	3 sec.	10	70						PI
2	90		15	70						
3	80		20	50						
4	90		25	60						
6	90		25	55	89	53	56	77	70	112
5	90		30	70						
7	80		30	70	30	15	28	15	25	14
8	80		30	70	55	40	41	30	25	13
9	80		35	70	35	25	40	25	43	26
10	80		40	45						
11	80		50	15						



Bull semen thawing ...

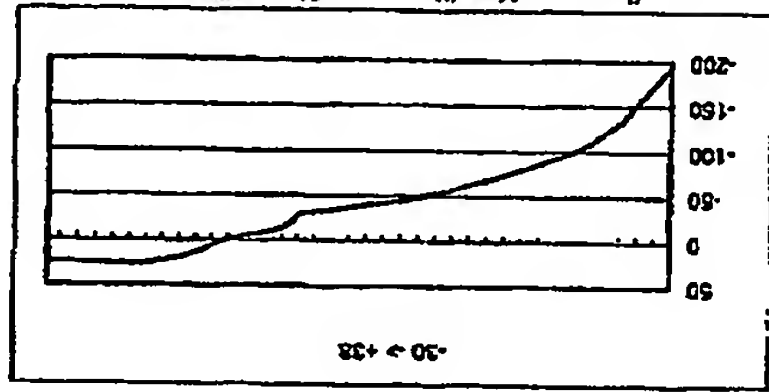


#	Algorithm	% motility
1	-186 > -24(20s) > 38	40
2	-186 > -24(20s) > 38	65
3	-186 > -24(20s) > 38	70
4	-186 > -24(20s) > 38	75
5	-186 > -24(40s) > 38	50
6	-186 > -24(40s) > 38	40
7	-186 > -24(40s) > 38	65
8	-186 > -24(40s) > 38	70

Average :

62.14286

1. Using water bath

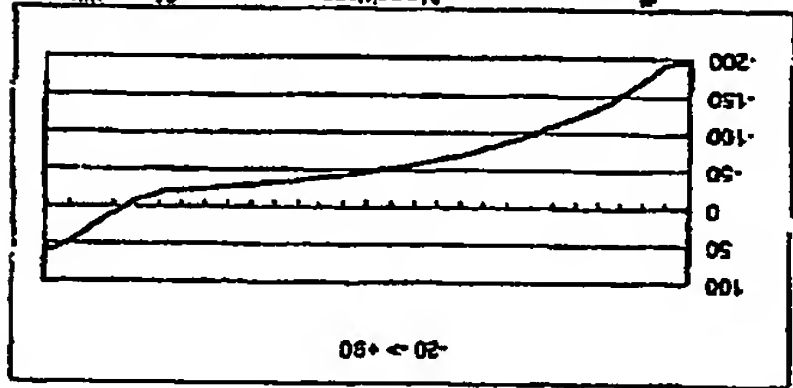

#	Algorithm	% motility
1	-186 > -24(20s) > 38	70
2	-186 > -24(20s) > 38	75
3	-186 > -24(20s) > 38	75
4	-186 > -24(20s) > 38	75
5	-186 > -24(40s) > 38	75
6	-186 > -24(40s) > 38	75
7	-186 > -24(40s) > 38	75

Average :

74.28571

2. Using new thawing machine

#	Alt time	Hot block	Temp.	Time (sec)	% motility
1	10	90	3	70	70
2	15	90	3	70	70
3	20	90	3	60	60
4	25	90	3	60	60
5	25	80	3	55	55
6	25	80	3	55	55
7	30	90	3	70	70
8	30	90	3	70	70
9	35	90	3	70	70

Stallion semen results

Stallion	Pre-freeze motility	Chilled semen after 30hrs motility	AO/PI %live	ORT	Planer Straw motility	AO/PI %live	ORT	MTG Tube motility	AO/PI %live	ORT
Nemrod	65	60	58.8	49.2	50	39.5	25.3	55	46.8	30.0
Jet set	80	70	59.4	47.1	65	54.0	41.6	65	76.5	50.5
William Curtis	50	50	62.3	36.0	20	25.4	18.1	40	55.0	37.7
Libra K	60	40	46.4	24.6	26	28.3	12.2	60	58.0	34.3
Samhite	60	70	70.0	59.1	60	51.4	33.0	70	68.0	40.3
Mean	57.0	58.0	67.3	43.2	42.0	39.5	26.0	58.0	61.1	38.6

Stallion semen results

Stallion	Pre-freeze motility	Chilled semen after 30hrs %motility AO/PI %live	ORT	Planer Straw %motility AO/PI %live	ORT	MTG Tube %motility AO/PI %live	ORT
Jat Set	70	46.4	24.9	50	49.3	37.1	60
Oboron	70	48.1	46.2	40	38.2	21.0	50
Rubek	80	68.2	57.9	60	54.0	24.0	60
Nemrod	60	64.1	52.1	30	34.3	12.0	50
Mean	70.0	60.0	45.3	45.0	44.2	23.5	55.0
		66.7					60.1
							31.7

Dear Udi

I hope this information is not too late, I have just received your e-mail. Please find attached the results for the latest MTG protocol for stallions. I am applying three tests for post thaw evaluation, namely Osmotic Resistance Testing (a membrane strength stress test), Acridine Orange/Propidium Iodide (a membrane viability test) and motility. None of these tests are unique to us. The freezing extender I am using contains the following:

Glucose monohydrate 15g  
Tri-sodium citrate 0.925g  
EDTA 0.925g  
Sodium hydrogen carbonate 0.3g  
Lactose 55g  
Lauryl sulphate 0.375g  
Lincospectin 1.0g  
Gentamycin 1.25ml  
Clarified egg yolk (centrifuged at 10000XG to remove fat) 200ml  
Glycerol 30ml (3%)  
Water to 1000ml

The clarified egg and low glycerol concentration makes this extender unique to us and could be regarded as specific or the MTG.

The freezing protocol I use is:

Manual seeding  
5°C start temp  
-50°C end temp  
1.0mm/second velocity.

A range of velocities can be applied from 0.6mm/second-3.0mm/second. These should all be protected. Also the start temp can be altered to a range between 25°C-5°C and the end temp can be between -5°C--100°C.

I think you are familiar with all other aspects (tubes etc.). Please do not hesitate to ring me if you need any other assistance. I am here for most of the day.

Regards  
Matt

The following section of this message contains a file attachment prepared for transmission using the Internet MIMB message format. If you are using Pegasus Mail, or any other MIME-compliant system, you should be able to save it or view it from within your mailer. If you cannot, please ask your system administrator for assistance.

Stallion	Pre-freeze motility	Chilled semen after 30hrs			Planer Straw			MTG Tube			Status
		%motility	AO/PI %live	ORT	%motility	AO/PI %live	ORT	%motility	AO/PI %live	ORT	
C. R. Gold	25	0	4	0	5	1	0	15	21.5	3.4	fail-fail
Libra-K	80	20	19.4	12.1	20	40.0	33.3	35	43.7	29.4	fail-pass
Samhite	80	60	72.4	46.2	3	31.9	23.3	60	48.6	48.0	fail-pass
Libra-K	80	60	79.1	48.4	20	42.8	33.9	40	57.0	46.0	fail-pass
Mill Law	40	30	42.3	39.2	10	34.5	28.2	30	42.2	35.8	fail-pass
Jester	90	60	71	59.1	30	29.3	20.1	40	40.3	23.2	fail-pass
Rob Roy	90	80	90.2	79.9	80	75.2	68.3	80	78.1	64.1	fail-pass
Pail Mail	70	20	18.4	0	25	32	16.3	35	38.4	26.2	fail-pass
Jester	60	50	48.6	41.2	30	34.2	22.6	50	53	43.4	fail-pass
Dramiro	70	50	55.2	40	20	24.3	11	40	49.4	41.5	fail-pass
Rubek	60	50	56.8	47.1	35	39.4	21.2	50	44.4	43.9	pass-pass
Rubek	60	60	52.3	51.1	50	49.8	35.9	60	49.2	40.2	pass-pass
Secundus	70	25	39.1	30.2	35	23.5	22	50	53.9	42.6	pass-pass
Schiller	90	60	71	49.8	60.0	43.5	34.7	60.0	64.5	51.2	pass-pass
Ludwig	80	40	39.4	28.4	35	35	28.9	40	36.8	29.4	pass-pass
Schiller	80	60	49.4	37.4	40	37.5	27.5	45	46.3	27.5	pass-pass
Secundus	80	50	62.8	50.1	40	43.2	36.2	50	51.4	37.9	pass-pass
C. R. Gold	0										
Mean	65.8	43.1	48.5	36.7	29.9	34.3	25.0	43.3	45.5	35.1	
Mean (fail-pass)	71.1	47.8	55.3	40.7	26.4	38.2	27.0	45.5	50.1	39.5	
Mean (pass-pass)	74.3	49.3	53.0	42.0	42.1	38.8	29.5	50.7	49.6	39.0	
					Straw	7/17 pass		Tube	16/17 pass		

27/06 '02 18:32 9872 3 562554

DR. W. FRIEDMAN ++ BTL

2005/011



Stallion semen results

To date 19.04.02


Stallion	Pre-freeze motility	Chilled semen after 30hrs			Pinger Straw			MTG Tube		
		%motility	AC/P1 %live	ORT	%motility	AC/P1 %live	ORT	%motility	AC/P1 %live	ORT
Samhine					5	31.9	21.3	60	48.6	48.0
Craigmanco Gold	60	0	4	0	5	1	0	15	21.6	3.4
Mill Law					10	34.5	26.2	30	42.2	35.8
Oberon	50	50	54.2		10	7.2	0	10	5.1	0
Beamus	60				10	26.6	12.6	45	41.1	35.1
Dramfro	70	50	58.2	40	20	24.3	11	40	49.4	41.5
Libra K	60				20	40.0	33.3	30	43.7	29.4
Libra K	60				20	42.6	33.6	40	37.0	46.0
William Curtis	50	50	52.3	38.0	20	23.4	18.1	40	55.0	37.7
Eagle	70				25	33	25	60	53.4	28.6
Eagle	70	60	68.4		25	18.1	12.1	60	52.6	49.6
Libra K	70				25	24.1	19	50	49.5	30.3
Libra K	60	40	48.4	24.6	25	26.3	12.2	60	58.0	34.3
Memphis	60				25	29.7	20.7	50	48.1	35.3
Nemrod	55	50	56.1		25	29.6	20.0	40	41.8	29.3
Pall Mall	70	20	15.4	0	25	32	16.3	35	38.4	28.2
Craigmanco Gold	70				30	28.4	27.8	50	44.6	29.8
Jester	80	60	71	58.1	30	28.3	20.1	40	40.3	23.2
Jester	60	50	48.6	41.2	30	34.2	22.6	50	53	43.4
Nemrod	60	60	64.1	52.1	30	34.3	12.0	50	49.7	21.6
Rubek	60	50	45.7		30	29.1	13.4	60	54.8	35.4
Rubek	60	50	53.5		30	34.6	15.4	60	56.8	41.1
Russel	70				30	21.9	17.4	80	57.1	39.4
Ludwig	80	40	39.4	28.4	35	35	28.8	40	36.6	29.4
Rob Roy	80	80	81.1	64	35	47.2	39.7	70	63.6	51.1
Rob Roy	90	80	71.8	62.7	35	41.4	32.2	50	58.2	34.6
Rubek	60				35	39.4	21.2	60	44.4	43.9
Secundus	70	25	38.1	30.2	35	23.6	22	5	53.9	61.2
Nemrod	70				40	44.0	28.6	60	50	30.7
Oberon	70	60	48.1	45.2	40	39.2	21.0	50	41.0	26.9
Schiller	60	50	49.4	37.4	40	37.5	27.5	45	48.3	27.5
Schiller	80	70	64.2	60	40	49.4	37.6	70	72.3	53.4
Schiller	80	70	69	58.4	40	49.7	35	60	68.4	47.3
Secundus	60	50	62.8	50.1	40	43.2	35.2	50	61.4	37.6
William	60	55	55.0		45	52.3	39.4	60	58.7	39.8
Jet Set	70	40	48.4	24.0	50	49.3	37.1	60	58.1	48.2
Memphis	60	50	47.9		50	53	29.6	60	57.4	37.1
Nemrod	65	60	58.5	49.2	50	38.5	25.3	55	46.6	30.0
Rubek	80				50	49.8	35.9	60	49.2	40.2
Samhine	80	70	70.0	59.1	50	51.4	33.0	70	68.0	40.3
William Curtis	70				50	49.7	23.0	60	58.4	44.1
Lagos	70				55	43	39.1	65	53.6	34.8
Lagos	60	60	68.1		55	57.6	28.6	60	63.4	50.0
Rubek	80	80	68.2	57.0	60	54.0	24.0	60	51.4	32.1
Schiller	90	60	71	49.8	60	43.5	34.7	60	64.5	51.2
Jet set	80	70	59.4	47.1	65	54.0	41.5	65	76.5	50.5
Rob Roy	90	80	80.2	79.0	80	73.2	56.3	80	76.1	64.1

27/08 '02 16:32 2872 3 682564 DR. M. FRIEDMAN \*\*\* BILL 2008/011

\_\_\_\_\_

Stallion	Pre-freeze mortality	Chilled semen after 30 hrs AOPI %live	Planor Straw AOPI %live	MTG Tube AOPI %live
Memphis	60	50	60	60
Oberon	50	50	53	57.4
Nemrod	55	50	25	41.5
Rubek 1	60	50	30	54.8
Rubek 2	60	50	30	60
Eagle	70	60	25	60
Lagos	60	60	30	60
William	60	55	45	60
Mean	59.4	53.1	33.8	48.8

temperature measured every second in the outer tube

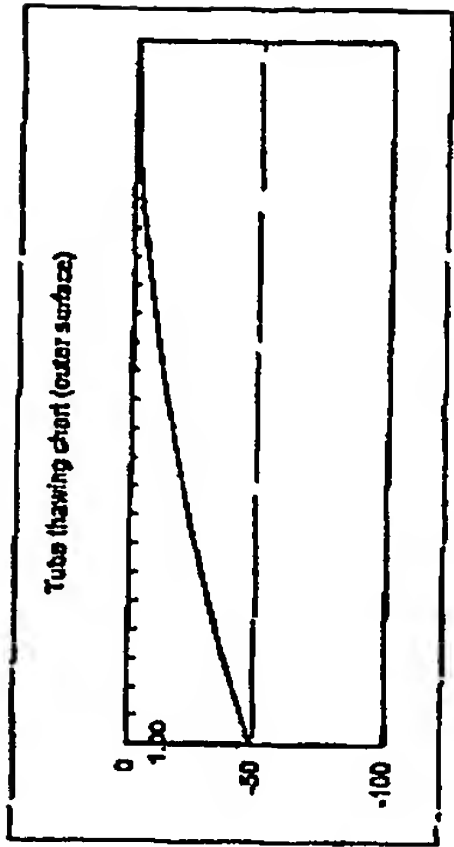

-48.5  
-49.3  
-49.1  
-48.9  
-48.7  
-48.5  
-48.3  
-48.1  
-47.9  
-47.7  
-47.5  
-47.3  
-47.1  
-46.9  
-46.7  
-46.5  
-46.3  
-46.1  
-45.9  
-45.7  
-45.5  
-45.3  
-45.1  
-45  
-44.9  
-44.7  
-44.6  
-44.4  
-44.1  
-44  
-43.8  
-43.6  
-43.4  
-43.2  
-43.1  
-42.8  
-42.6  
-42.4  
-42.3  
-42.1  
-41.9  
-41.7  
-41.5  
-41.4  
-41.2  
-41.1  
-40.9  
-40.8  
-40.6  
-40.4  
-40.2  
-40.1  
-39.8  
-39.7



-39.5  
-39.3  
-39.2  
-39  
-38.8  
-38.7  
-38.5  
-38.4  
-38.2  
-38.1  
-37.9  
-37.7  
-37.6  
-37.4  
-37.4  
-37.2  
-37  
-36.9  
-36.8  
-36.8  
-36.4  
-36.3  
-36.2  
-36  
-35.8  
-35.7  
-35.6  
-35.4  
-35.3  
-35.1  
-34.9  
-34.8  
-34.7  
-34.5  
-34.4  
-34.2  
-34.1  
-33.9  
-33.8  
-33.7  
-33.6  
-33.5  
-33.3  
-33.1  
-33  
-32.9  
-32.7  
-32.6  
-32.4  
-32.3  
-32.1  
-32  
-31.9  
-31.7  
-31.5  
-31.4



-31.3  
-31.1  
-31  
-30.8  
-30.7  
-30.5  
-30.4  
-30.2  
-30.1  
-30  
-29.8  
-29.7  
-29.6  
-29.5  
-29.3  
-29.1  
-29  
-28.9  
-28.7  
-28.6  
-28.4  
-28.3  
-28.2  
-28  
-27.8  
-27.8  
-27.5  
-27.5  
-27.3  
-27.2  
-27.1  
-26.9  
-26.7  
-26.8  
-26.5  
-26.4  
-26.3  
-26.1  
-26  
-25.8  
-25.6  
-25.5  
-25.3  
-25.2  
-25  
-24.9  
-24.8  
-24.6  
-24.5  
-24.3  
-24.2  
-24.1  
-24  
-23.8  
-23.7  
-23.6



-23.5  
-23.3  
-23.2  
-23  
-22.9  
-22.8  
-22.8  
-22.6  
-22.5  
-22.3  
-22.2  
-22.1  
-22  
-21.9  
-21.7  
-21.6  
-21.5  
-21.4  
-21.3  
-21.1  
-21.1  
-20.9  
-20.8  
-20.7  
-20.5  
-20.4  
-20.3  
-20.2  
-20.1  
-20  
-19.9  
-19.8  
-19.6  
-19.5  
-19.4  
-19.4  
-19.2  
-19.2  
-19.1  
-18.9  
-18.8  
-18.7  
-18.6  
-18.5  
-18.4  
-18.2  
-18.1  
-18  
-17.9  
-17.8  
-17.7  
-17.6  
-17.5  
-17.4  
-17.2  
-17.1

-17  
-16.9  
-16.8  
-16.7  
-16.6  
-16.5  
-16.4  
-16.3  
-16.2  
-16  
-15.9  
-15.8  
-15.7  
-15.6  
-15.5  
-15.4  
-15.3  
-15.2  
-15.1  
-15  
-14.8  
-14.7  
-14.6  
-14.5  
-14.4  
-14.3  
-14.2  
-14.1  
-14  
-13.9  
-13.8  
-13.7  
-13.6  
-13.5  
-13.4  
-13.3  
-13.2  
-13.1  
-13  
-12.9  
-12.8  
-12.7  
-12.6  
-12.5  
-12.4  
-12.3  
-12.2  
-12.2  
-12.1  
-12  
-11.9  
-11.8  
-11.8  
-11.7  
-11.6



-11.5  
-11.4  
-11.3  
-11.2  
-11.1  
-11.1  
-11  
-10.9  
-10.7  
-10.7  
-10.6  
-10.5  
-10.5  
-10.4  
-10.3  
-10.1  
-10  
-10  
-9.9  
-9.9  
-9.7  
-9.7  
-9.6  
-9.6  
-9.4  
-9.3  
-9.2  
-9.2  
-9  
-8.9  
-8.8  
-8.8  
-8.7  
-8.7  
-8.6  
-8.5  
-8.4  
-8.3  
-8.2  
-8.1  
-8  
-8  
-7.8  
-7.7  
-7.7  
-7.6  
-7.5  
-7.4  
-7.3  
-7.2  
-7.2  
-7.1  
-7  
-6.9  
-6.9  
-6.8  
-6.7



-8.8  
-8.6  
-8.5  
-8.4  
-8.3  
-8.2  
-8.1  
-8  
-8  
-5.9  
-5.8  
-5.7  
-5.7  
-5.6  
-5.5  
-5.4  
-5.4  
-5.3  
-5.3  
-5.2  
-5.1  
-5  
-4.9  
-4.9  
-4.8  
-4.8  
-4.7  
-4.6  
-4.5  
-4.5  
-4.4  
-4.3  
-4.3  
-4.2  
-4.1  
-4  
-4  
-3.9  
-3.8  
-3.7  
-3.7  
-3.6  
-3.6  
-3.5  
-3.4  
-3.4  
-3.3  
-3.2  
-3.1  
-3.1  
-3  
-2.9  
-2.9  
-2.8  
-2.8

-2.8  
-2.7  
-2.7  
-2.6  
-2.6  
-2.5  
-2.5  
-2.4  
-2.4  
-2.3  
-2.3  
-2.2  
-2.1  
-2.1  
-2.1  
-2  
-2  
-2  
-1.9  
-2  
-1.8  
-1.9  
-1.9  
-1.9  
-1.9  
-1.9  
-1.8  
-1.8  
-1.8  
-1.8  
-1.7  
-1.7  
-1.6  
-1.6  
-1.6  
-1.6  
-1.5  
-1.5  
-1.5  
-1.5  
-1.5  
-1.4  
-1.4  
-1.4  
-1.4  
-1.4  
-1.4  
-1.4  
-1.3  
-1.3  
-1.3  
-1.3  
-1.3  
-1.3  
-1.3






[illegible]

temperature measured every second in the inner tube

-49.5  
-49.2  
-49  
-48.8  
-48.6  
-48.4  
-48.2  
-48  
-47.8  
-47.6  
-47.4  
-47.2  
-47  
-46.7  
-46.5  
-46.2  
-46  
-45.8  
-45.6  
-45.4  
-45.2  
-45  
-44.8  
-44.5  
-44.4  
-44.2  
-44  
-43.7  
-43.6  
-43.4  
-43.2  
-43  
-42.8  
-42.6  
-42.4  
-42.2  
-42  
-41.8  
-41.6  
-41.5  
-41.3  
-41  
-40.9  
-40.7  
-40.6  
-40.4  
-40.2  
-40.1  
-39.8  
-39.7  
-39.5  
-39.3  
-39.2  
-39

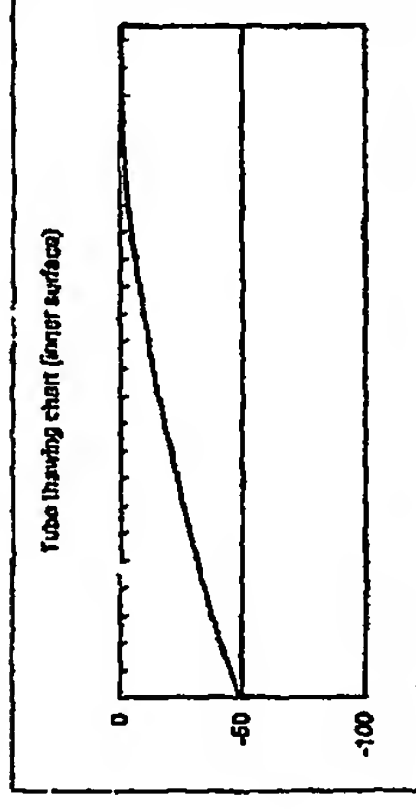

-38.9  
-38.7  
-38.5  
-38.3  
-38.2  
-38  
-37.8  
-37.6  
-37.5  
-37.3  
-37.1  
-37  
-36.9  
-36.7  
-36.6  
-36.4  
-36.3  
-36.1  
-35.9  
-35.8  
-35.6  
-35.5  
-35.3  
-35.1  
-34.9  
-34.8  
-34.6  
-34.5  
-34.3  
-34.1  
-34  
-33.8  
-33.6  
-33.4  
-33.3  
-33.2  
-33.2  
-33.1  
-32.9  
-32.8  
-32.7  
-32.5  
-32.4  
-32.3  
-32.1  
-31.9  
-31.8  
-31.7  
-31.6  
-31.4  
-31.2  
-31  
-30.9  
-30.7  
-30.5  
-30.4

-30.2  
-30.1  
-29.9  
-29.7  
-29.6  
-29.6  
-29.4  
-29.3  
-29.1  
-28.9  
-28.8  
-28.6  
-28.6  
-28.4  
-28.2  
-28.1  
-28  
-27.9  
-27.7  
-27.5  
-27.4  
-27.2  
-27  
-27  
-26.8  
-26.6  
-26.5  
-26.4  
-26.1  
-26  
-25.9  
-25.8  
-25.7  
-25.5  
-25.4  
-25.2  
-25.1  
-25  
-24.8  
-24.7  
-24.6  
-24.4  
-24.3  
-24.2  
-24  
-23.9  
-23.8  
-23.6  
-23.5  
-23.4  
-23.2  
-23.2  
-23  
-22.9  
-22.7  
-22.6

-22.4  
-22.4  
-22.3  
-22.1  
-22  
-21.9  
-21.7  
-21.6  
-21.5  
-21.4  
-21.3  
-21.1  
-21  
-20.9  
-20.8  
-20.7  
-20.5  
-20.4  
-20.3  
-20.2  
-20.1  
-19.9  
-19.8  
-19.7  
-19.5  
-19.4  
-19.3  
-19.2  
-19.1  
-18.9  
-18.9  
-18.8  
-18.6  
-18.5  
-18.4  
-18.3  
-18.2  
-18.1  
-17.9  
-17.8  
-17.7  
-17.6  
-17.5  
-17.4  
-17.3  
-17.2  
-17.1  
-16.9  
-16.8  
-16.7  
-16.6  
-16.4  
-16.3  
-16.2  
-16.1  
-16



-15.8  
-15.7  
-15.5  
-15.4  
-15.4  
-16.3  
-15.1  
-16  
-14.9  
-14.8  
-14.6  
-14.5  
-14.4  
-14.3  
-14.2  
-14.1  
-13.9  
-13.8  
-13.7  
-13.6  
-13.5  
-13.4  
-13.3  
-13.2  
-13.1  
-13  
-12.9  
-12.8  
-12.7  
-12.6  
-12.4  
-12.4  
-12.3  
-12.2  
-12.1  
-12  
-11.9  
-11.9  
-11.8  
-11.7  
-11.6  
-11.5  
-11.4  
-11.3  
-11.2  
-11.1  
-11  
-10.9  
-10.8  
-10.7  
-10.6  
-10.5  
-10.4  
-10.3  
-10.3  
-10.2



-10.1  
-9.9  
-9.8  
-9.7  
-9.6  
-9.5  
-9.4  
-9.3  
-9.2  
-9.1  
-9  
-8.9  
-8.9  
-8.8  
-8.7  
-8.6  
-8.4  
-8.4  
-8.3  
-8.2  
-8.1  
-8.1  
-8  
-7.9  
-7.8  
-7.7  
-7.6  
-7.5  
-7.4  
-7.3  
-7.2  
-7.1  
-7  
-6.9  
-6.9  
-6.7  
-6.7  
-6.6  
-6.5  
-6.4  
-6.3  
-6.3  
-6.2  
-6.1  
-6  
-5.9  
-5.9  
-5.8  
-5.7  
-5.6  
-5.5  
-5.5  
-5.4  
-5.3  
-5.2  
-5.1

-5.1  
-5.1  
-5  
-4.9  
-4.8  
-4.7  
-4.7  
-4.6  
-4.6  
-4.5  
-4.4  
-4.3  
-4.3  
-4.1  
-4  
-4  
-3.9  
-3.8  
-3.8  
-3.7  
-3.6  
-3.6  
-3.5  
-3.4  
-3.4  
-3.3  
-3.2  
-3.2  
-3.1  
-3.1  
-3  
-2.9  
-2.8  
-2.8  
-2.8  
-2.7  
-2.7  
-2.6  
-2.5  
-2.5  
-2.5  
-2.5  
-2.4  
-2.4  
-2.3  
-2.2  
-2.2  
-2.2  
-2.1  
-2.1  
-2.1  
-2.1  
-2  
-2  
-1.9  
-1.9

[illegible]


0.9 0.9

[illegible]

UNITED STATES PATENT AND TRADEMARK OFFICE  
DOCUMENT CLASSIFICATION BARCODE SHEET



**Oath/Declaration, Small Entity,  
and Power of Attorney**

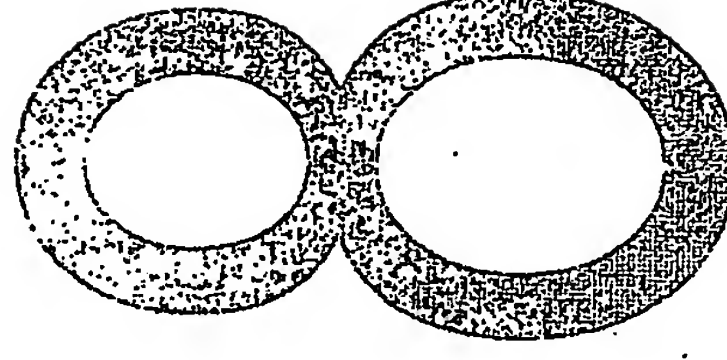

Level - 2  
Version 1.1  
Updated - 8/01/01



Form PTO-1585 (Rev. 5-92) OMB No. 051-0011 (exp. 4/94)		RECORDATION FORM COVER SHEET <b>PATENTS ONLY</b>		U.S. DEPARTMENT OF COMMERCE Patent and Trademark Office	
Tab settings: =>=>=>					
To the Honorable Commissioner of Patents and Trademarks: Please record the attached original document or copy thereof.					
1. Name of conveying party(ies): AMIR ALAY WEIR URI		2. Name and address of receiving party(ies): Name: <u>Interface Multigrid Technology Ltd.</u> Internal Address: _____ Street Address: <u>3 Hamazma St.</u> <u>Ness Ziona 70400, Israel</u> City: _____ State: _____ ZIP: _____ Additional name(s) & address(es) attached? <input type="checkbox"/> Yes <input checked="" type="checkbox"/> No			
Additional name(s) of conveying party(ies) attached? <input type="checkbox"/> Yes <input checked="" type="checkbox"/> No					
3. Nature of conveyance: <input checked="" type="checkbox"/> Assignment <input type="checkbox"/> Merger <input type="checkbox"/> Security Agreement <input type="checkbox"/> Change of Name <input type="checkbox"/> Other: _____ Execution Date: <u>27 June 2002</u>					
4. Application number(s) or patent number(s): If this document is being filed together with a new application, the execution date of the application is: <u>27-Jun-02</u> A. Patent Application No.(s) _____ B. Patent No.(s) _____ Additional numbers attached? <input type="checkbox"/> Yes <input checked="" type="checkbox"/> No					
5. Name and address of party to whom correspondence concerning document should be mailed: Name: <u>Mark M. Friedman</u> Internal Address: _____ Street Address: <u>c/o Bill Reikenghorn</u> <u>Discovery Biotech, Florin Way 9003</u> City: <u>Upper Merion, State: MD ZIP: 20773</u>		6. Total Number of applications and patents involved: <u>1</u>			
		7. Total fee (37 CFR 3.41).....\$ <u>40</u> <input type="checkbox"/> Enclosed <input checked="" type="checkbox"/> Authorized to be charged to deposit account			
		8. Deposit account number: <u>06-2140</u>			
		(Attach duplicate copy of this page if paying by deposit account)			
		DO NOT USE THIS SPACE			
9. Statement and signature. To the best of my knowledge and belief, the foregoing information is true and correct and any attached copy is a true copy of the original document. Mark M. Friedman Name of Person Signing		Signature <u>[Signature]</u> Date <u>27 JUN 02</u> Total number of pages including cover sheet, attachments, and document <u>2</u>			

Mail documents to be recorded with required cover sheet information to:  
Commissioner of Patents & Trademarks, Box Assignments  
Washington, D.C. 20231

WPS

**ASSIGNMENT**

For good and valuable consideration, the receipt and sufficiency of which is hereby acknowledged, the undersigned:

Amir Arav                      Meir Uri

(hereinafter called the "assignor(s)"), hereby sell(s), assign(s) and transfer(s) to:

Interface Multigrad Technology Ltd.  
3 Hamazmera St.  
Ness Ziona 70400  
Israel

(hereinafter called the "assignee(s)"), its/their successors, assigns, nominees or other legal representatives, the Assignor's entire right, title and interest in and to the invention entitled:


**EMBRYO GUARD**

described and claimed in the following patent applications:

U.S. Provisional Application identified as Attorney docket No. 791/14 and executed the same date as this assignment;

and in and to said Patent Applications, and all original and reissued Patents granted therefor, and all divisions and continuations thereof, including the right to apply and obtain Patents in all other countries, the priority rights under International Conventions, and the Letters Patent which may be granted thereon;

Signed and sealed this 22 day of June 2002

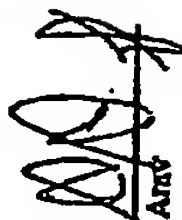  
Amir Arav

  
Meir Uri



16.05.02-19.05.02

Station	Pre-freeze molality	Chilled semen after 30hrs %molality AO/P1 %live	Planer Straw 7/17 pass %molality AO/P1 %live	MTO Tube 16/17 pass %molality AO/P1 %live	Status
C. R. Gold	25	0	5	15	tall-fail
Libra-K	80	19.4	20	29.4	fail-pass
Samhire	60	72.4	3	35	fail-pass
Libra-K	80	79.1	20	40	fail-pass
Mill Low	40	42.3	10	30	fail-pass
Jester	80	71	30	40	fail-pass
Rob Roy	80	80.2	80	23.2	fail-pass
Pall Mall	70	16.4	25	80	fail-pass
Jester	60	48.6	30	35	fail-pass
Dramiro	70	56.2	20	53	fail-pass
Rubek	60	58.8	35	40	tail-pass
Rubek	60	52.3	60	50	pass-pass
Secundus	70	62.3	50	60	pass-pass
Schiller	80	39.1	35	60	pass-pass
Ludwig	80	71	60.0	60.0	pass-pass
Schiller	80	49.4	36	64.5	pass-pass
Secundus	80	28.4	40	38.6	pass-pass
C. R. Gold	0	50.1	40	51.4	pass-pass
Mean	65.8	43.1	38.7	43.3	Mean
Mean (fall-pass)	71.1	47.8	26.4	45.6	Mean (fall-pass)
Mean (pass-pass)	74.3	49.3	42.0	50.7	Mean (pass-pass)
			Straw 7/17 pass	Tube 16/17 pass	
			42.1	50.7	
			38.8	49.5	
			38.2	50.1	
			27.0	39.5	
			34.3	48.5	
			25.0	35.1	

202081 2020 01

**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

☒ **BLACK BORDERS**

☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**

☐ **FADED TEXT OR DRAWING**

☐ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**

☐ **SKEWED/SLANTED IMAGES**

☒ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**

☐ **GRAY SCALE DOCUMENTS**

☐ **LINES OR MARKS ON ORIGINAL DOCUMENT**

☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**

☐ **OTHER:** \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**